



T-Max Hospital Autoclaves

For CSSDs, OR and Medical Centers



T-Max Line


Tuttnauer



Hospital Autoclaves

Tuttnauer hospital autoclaves ensure reliable sterile processing for Hospital CSSD, OR, and Medical Centers. The sterilizers are designed and manufactured in a state-of-the-art facility in compliance with strict international standards to ensure safe re-use of sterilized equipment.

Experience Since 1925


Tuttnauer has manufactured products for over 100 years that have developed a reputation for their quality, high performance and reliability, sophisticated features, and ability to satisfy customer expectations. Tuttnauer's sterilization & infection control products are trusted at over 350,000 installations worldwide including Hospitals, Clinics and Laboratories.

Exceptional Support Worldwide

Tuttnauer's staff provides expert pre-sales and after-sales support services to satisfy customer expectations. Our teams are made up of multi-cultural people who are able to comfortably work with customers around the world. Tuttnauer provides in-depth face-to-face training in many locations around the world to ensure that Tuttnauer technicians and engineers are experts in their ability to support each customer's technical service needs.

At Tuttnauer we highly value customer feedback which contributes to the continuous improvement of our products and support services.



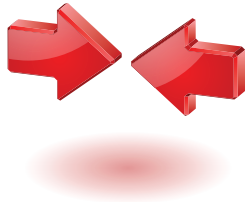

Narrow Design for Limited Space

A narrow design autoclave is as wide as 99cm and can sterilize a fully loaded chamber with 6 to 8 STU. Facilities with space layout limitations will benefit from narrow design autoclaves which provide high capacity in a limited space. Narrow design autoclaves allow for more autoclaves in a given space to ensure continued sterile processing in case of autoclave downtime.




Uniform Heat Distribution

Fully jacketed chambers ensure uniform heat distribution to ensure consistent sterilization conditions for all loads inside the chamber. A chamber that is not fully jacketed may have cool spots which can cause improper sterilization.



Stainless Steel Piping

Stainless steel piping and connectors are used to avoid corrosion which is necessary to safeguard the integrity of the autoclave and prevent contamination of the sterilized loads (either as chemical or particulate contamination).

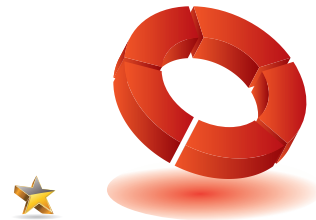

Perfectly Optimized Cycles

With years of research and customer feedback we have fine-tuned our formula for perfectly optimized cycles with the precise combination of pulses, timing, temperature control, speed and drying.



Load Traceability

The optional R.P.C.R software allows for tracing sterilized loads by barcode. All cycle records and associated barcodes are automatically recorded to a PC on your network. In addition, R.P.C.R allows for convenient access to cycle reports, including graphs and tables (in PDF format).

Advanced Control System for Your CSSD

Take advantage of Tuttnauer's sophisticated user-friendly control systems for repeatable high performance.

Standard Features

- 7" Multi-color touch screen panel
- Keypad control panel on second door of two door autoclaves with Bacsoft controller
- Stores the last 200 cycles in built-in memory (Bacsoft)
- Multiple access levels and user passwords to control access/operation of the autoclave
- Diagnostic In/Out test (enables technician to check each system component separately)
- Sterilization Temperature range 105°C to 138°C

Load Tracing with Barcodes

When loading the autoclave a barcode reader can be used to scan the barcodes of each load. Once the cycle starts the barcode numbers will be printed together with cycle information.

Optional Features

- 21 CFR part 11

Process Workflow Management Data Connection


Provides real-time data to any process workflow management system supporting Modbus

Cycle Programs

Horizontal Autoclaves are preprogrammed with 8 perfectly optimized cycle programs that cover the processing needs of most CSSDs and Operating Rooms.

8 Preset Cycle Programs

- Unwrapped cycle for operating rooms (OR) that need surgery instruments for immediate use
- Wrapped and double wrapped loads for CSSD processing on instruments and textiles at 134°C and 121°C for delicate instruments
- Additional special program cycles for prion applications and special wrapping materials

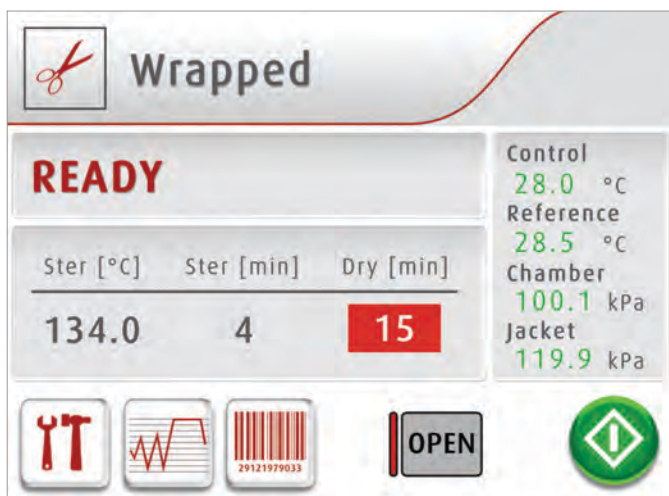

Sophisticated Touch Screen HMI

The HMI (Human Machine Interface) has been designed with the following considerations:

- Multi-color display for easier reading from a distance
- Multilingual (26 languages)
- Graphical display of Temperature and Pressure trend graphs




2 Test Cycles

- Bowie Dick/Helix for steam penetration testing
- Leak test to check vacuum integrity

Custom Cycles

- 20 cycles available for customized cycle programs

Bacsoft Controller

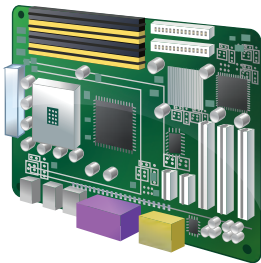

STANDARD

OPTIONAL




Printer



Touch Screen
Loading side



Remote Internet access
for Technical Support
(requires R.P.C.R. and local
cellular router for
Internet connection)

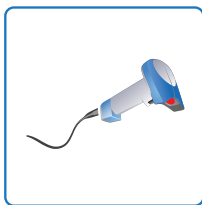

Barcode
Reader



R.P.C.R software
for access from
remote PC on
internal network

Load Tracing & Automatic Recording of Cycle Information to Your PC

R.P.C.R Software (optional)


Cycle Records You Can Rely On

- Automatic recording of cycle information to any PC on your network via Ethernet
- Convenient access to graphs and tables that are easy to understand
- Easily generate PDF reports
- No need to file printouts, saving you time

Load Traceability

Conveniently trace loads from a PC on your network with R.P.C.R software. When barcodes of loads are scanned before a cycle starts, the barcode information is reported together with the cycle records.

With R.P.C.R you can see: Graphs of the cycle data, Numeric cycle data, Cycle print-outs, Measured values table, Parameter table.

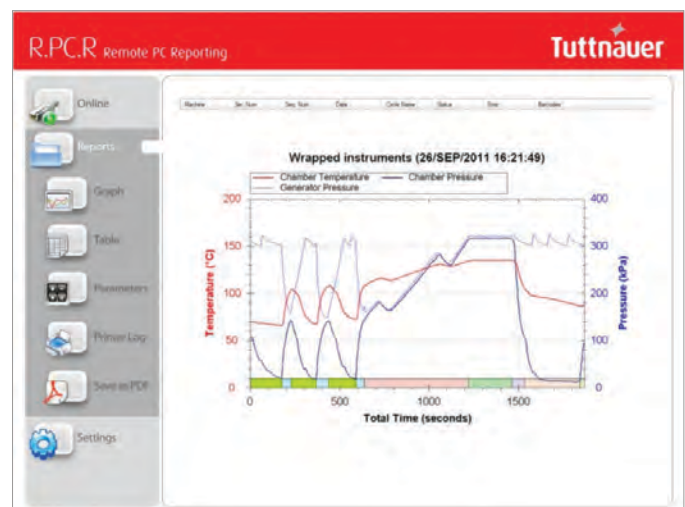

T-Max Hospital Autoclaves

The T-Max Narrow Series

6 to 10 StU Capacity

T-Max Narrow Series Hospital Autoclaves have chamber volumes ranging in size from 430 Liters to 700 Liters (6 to 10 StU) and are available with fully automatic vertical sliding doors.



Model	Chamber Volume (Liter)	StU Capacity	Chamber Dimensions (WxHxD) mm
T-Max 6	430 Liter	6	660 x 660 x 990
T-Max 8	565 Liter	8	660 x 660 x 1295
T-Max 10	700 Liter	10	660 x 660 x 1620

Model 1V/2V	EXTERNAL DIMENSIONS (mm)			
	Width	Height	Depth	
	1 & 2 Doors	1 & 2 Doors	1 Door	2 Doors
T-Max 6	999	1980	1400	1310
T-Max 8	999	1980	1690	1610

Model 2V	EXTERNAL DIMENSIONS (mm)		
	Width	Height	Depth
	2 Doors	2 Doors	2 Doors
T-Max 10	999	1980	2010



Vertical Sliding Door



Loading Equipment

High quality stainless steel loading equipment for loading and unloading.

- Pull Out Trays

Stainless steel trays equipped with rails for easy loading and unloading. The rails are designed to prevent the trays from rolling over.

- Loading Carts and Transfer Carriages

The 316L loading carts are designed to roll from the transfer carriage onto the chamber rails for easy handling of heavy loads. To ensure safety and ease-of-use the carriage is equipped with a lock preventing sliding of the cart. Swivel wheels with wheel breaks maximize mobility in limited space.



Complimentary Products

Ultrasonic Cleaner

Deep cleaning of hollows before sterilization



Sealing Machine

For sealing instruments in pouches for sterile storage



Safety

Safety for personnel, autoclave and load are priority in the design, construction and operation of any Tuttnauer autoclave. Tuttnauer is committed to the highest industry safety standards and directives to ensure safety not only for your employees operating the autoclaves but also for your sterile processing facility and the loads being sterilized.

Tuttnauer autoclaves are provided with redundant independent monitoring systems and audio visual alarms to notify operators of any issue that requires attention. An emergency stop button on the loading side of the autoclave may be used to safely stop the autoclave cycle.

Company standards and legislation

Directives, Guidelines & Technical Standards

- ISO 13485:2016 Quality Management System – Medical devices
- Medical Device Regulation 2017/745
- Medical Device Single Audit Program- companion document (MDSAP)
- ISO 9001: 2015 “Quality Management System”
- ISO 14001:2015 environmental management system
- 2012/19/EU WEEE and EU 2017/2102 Restriction of the use of certain hazardous substances (RoHS Directive)

Product Standards

Pressure Vessel- ASME Code section I and section VIII. Div.I
Pressure Vessel- PED 2014/68/EU
Pressure Vessel- Chinese Regulations- Special Equipment Licensing Office
EN 285 “Large steam sterilizer”

ANSI/AAMI/ST8 Hospital Steam Sterilizers

IEC 61010-1 Safety requirements for electrical equipment for measurement, control, and laboratory use - Part 1: General requirements

EN IEC 61010-2-040: Safety requirements for electrical equipment for measurement, control, and laboratory use - Part 2-040: Requirements for sterilizers and washer-disinfectors used to treat medical materials

EN 61326-1 EMC Requirements for Electrical Equipment
IEC 62304: Medical Device Software; Software life cycle processes

Tuttnauer pressure vessels are both ASME and PED certified. All ASME certified vessels are inspected by an independent authorized ASME inspector.



Learn from our Experts

Join our blog tuttnauer.com/blog



Hospital Sterilizers

International Sales and Marketing
Contact us at <https://tuttnauer.com/contact>
Website: tuttnauer.com



Tuttnauer Europe b.v.
Hoeksteen 11, 4815 PR
P.O.B. 7191, 4800 GD Breda
The Netherlands
Tel: +31 (0) 765 423 510 Fax: +31 (0) 765 423 540
E-mail: info@tuttnauer.nl

Tuttnauer USA Co.
345 Oser Ave
Hauppauge, NY11788
Tel: +800 624 5836, +631 737 4850 Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com



Distributed by: