


Giving

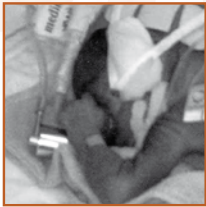

future...

Miniflow nCPAP System

- nasal CPAP-generator,
compatible with ventilators
- noiseless



Giving future...



medin
Medical Innovations

Functions & advantages

Miniflow nCPAP

for single use

- A) Inspiration
- B) Expiration
- C) Separation wall
- D) Flexible prong adapter

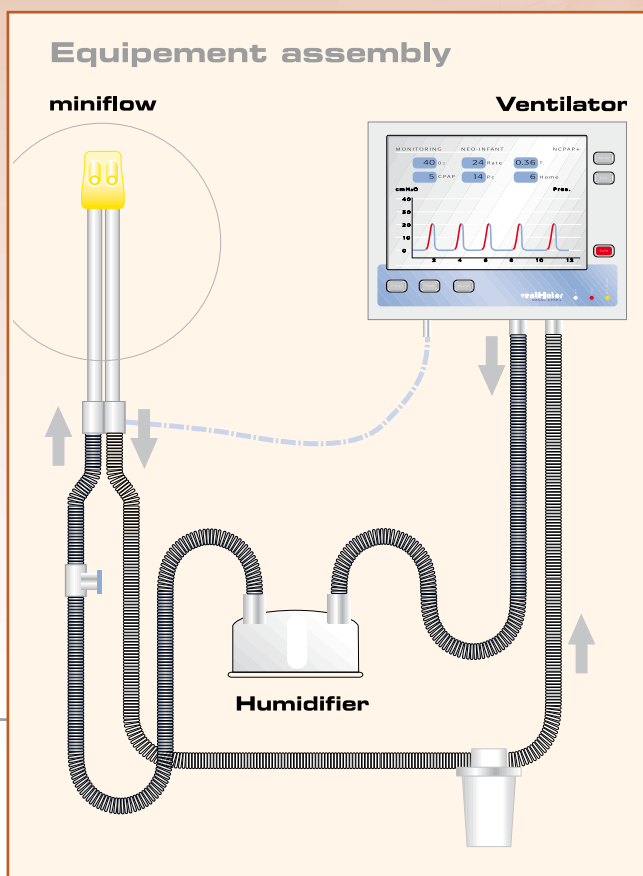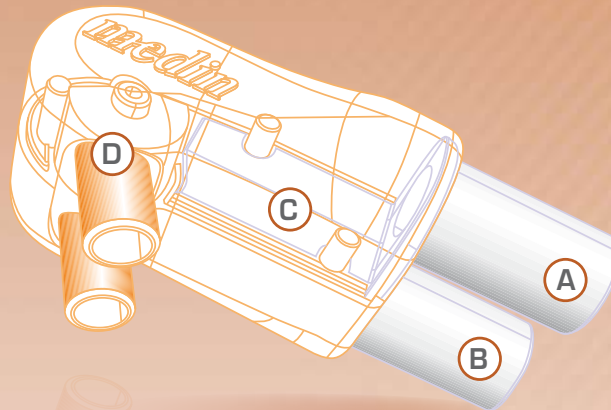

Advantages of the Miniflow nCPAP system:

- Minimal dead space
- Separation of inspiration- and expiration-gases
- Compatible with Medijet® accessoires
- Reduction of work of breathing
- Flexible prong adapter for best position
- Noiseless
- Simple use



- Always use the largest prong to avoid leakage
- Switch between prong and mask after some hours
- Insert the prong only halfway into the nose

IMPORTANT: Single use product for a single patient only!

Miniflow is a patient interface for the CPAP mode of ventilators. Miniflow can therefore be connected to the inspiration and expiration tubes and thus to the ventilator. The ventilator adjusts and controls the CPAP parameter and provides the necessary gas flow.

Connect Miniflow to the inspiration and the expiration tube, if your ventilator needs an external pressure measurement, add the T-piece into this connection (REF 4010). Close the prong adapter of the Miniflow with your finger and adjust the nCPAP parameters on the ventilator. If this is not possible check the connections to see if one is leaking.

Choose one prong or mask to fit the to patient's nose, connect the selected prong to Miniflow, the loops must point to the opposite direction of the tubes.

Choose the right bonnet, the measurement strip (REF 2150) shows you the necessary size.

Open the hook and loop fastener, slide the strips through the two holes in the prong and mask

Choose the right angle on Miniflow, insert the prong or mask with the adjusted Miniflow into the baby's nostrils, fix the strips with the small hooks to the bonnet below the ears and never over the eyes.

Take care that the strips are not tightened too much. Position Miniflow's tubes in the butterfly and close the butterfly. Afterwards check the position of the bonnet and strips again.

See user manual of Miniflow.

Miniflow nasal CPAP System for neonates & premature infants



Other circuits available

Miniflow & Circuits:

Complete Set	Part No./ REF	Part No./REF
10 x Miniflow & circuits (without chamber) MR850	40206749	40206748
Miniflow alone		
20 pieces	4000	
T-piece pressure tube	4010	
Circuits alone		
10 sets, without chamber	206749 heated exp.	206748 water trap in exp.
10 sets, chamber included	206747 heated exp.	206746 water trap in exp.



Masks & Prongs:

Starter Kit Masks & Prongs	Part No./ REF
20 masks and prongs in various sizes	1000-50

Prongs

Size	Part No./ REF	Size	Part No./ REF
micro	1200-01	large	1200-22
small	1200-21	x-large	1200-03
medium	1200-02	Asia ¹	1200-33

Masks

Size	Part No./ REF	Size	Part No./ REF
small	1200-04	large	1200-06
medium	1200-05		

¹ wide version of L



Bonnets:

Starter Kit Single Use Bonnets	Part No./ REF
14 bonnets, 2 of every size	1000-60

single use

Size	Part No./ REF	Size	Part No./ REF
XS	1214-10	XL	1218-10
S	1215-10	XXL	1219-10
M	1216-10	XXXL	1220-10
L	1217-10		

re-use

Size	Part No./ REF	Size	Part No./ REF
XS	1214	XL	1218
S	1215	XXL	1219
M	1216	XXXL	1220
L	1217	XXXXL	1221

Measuring tape:

Part No./ REF
2150



Distributors worldwide



Your distributor:

Austria · Australia · Azerbaijan · Bangladesh · Belarus
Belgium · Canada · China
Combined Arab Emirates
Egypt · Czech Republic
Denmark · Estonia · Finland
France · Germany · Greece
Great Britain · Hong Kong
Hungary · India · Indonesia
Iran · Israel · Italia · Japan
Jordan · Kazakhstan
Kirgistan · Korea · Kuwait
Latvia · Lithuania · Luxembourg · Malaysia · Morocco
Mexico · Moldavia · Netherlands · Norway · Pakistan
Peru · Poland · Rumania
Russia · Saudi Arabia
Singapore · Slovenia · Slovakia · South Africa · Spain
Sri Lanka · Sweden
Switzerland · Syria · Thailand
Taschikistan · Turkey · Ukraine

www.medin-innovations.com

Medical Innovations GmbH
Lindberghstr. 1
D-82178 Puchheim

Phone: +49 89 89670950
Fax: +49 89 80048476
info@medin-innovations.com

medin
Medical Innovations



Certification in accordance with DIN EN ISO 9001 and DIN EN ISO 13485,
Annex II section 3 of the Council Directive 93/42/EEC