



SKB039(H)

STRETCHER TROLLEY

Registration code: SKC-3B5



SAIKANG Medical

Healthier
your life

SKB039(H)

FEATURES AND HIGHLIGHTS



Technical parameters

▪ Unfoldingsize(L×W×H):	1950×570×910mm(±10mm)
▪ Foldingsize(L×W×H):	2030×570×300mm(±10mm)
▪ Heightadjustment:	78cm,87cm,91cm
▪ Backrestangle:	0-75°(±5°)
▪ Legangle:	0-15°(±5°)
▪ H-shaped legs	
▪ Load capacity	≤ 220KG
▪ One staff member is required to operate	
▪ The stretcher surface can be separated and has a locking device.	
▪ Aluminum alloy, stainless steel	

Technical configuration

☑ Casters:	6-inch wheels × 4 (2 brakes)
☑ Strap:	1 chest strap; One leg strap;
☑ Mattress(Thickened sponge waterproof PVC, thickness 100mm):	1pc
☑ Scalable saltwater rack:	1 unit
☑ 2-stage design, with a maximum adjustable length of 71cm	
☑ Stretcher fixing lock:	1 set

Design

In line with ergonomics, it meets the transfer requirements of patients, and its design emphasizes usability and multifunctionality.




Steel

It is made of high-strength aluminum alloy material, characterized by its compact size, light weight, durability, and resistance to rust. It has a maximum weight-bearing capacity of up to 180 kilograms.




Welding

The argon gas shielding welding process ensures the creation of well-formed weld seams and a smooth, clean surface.

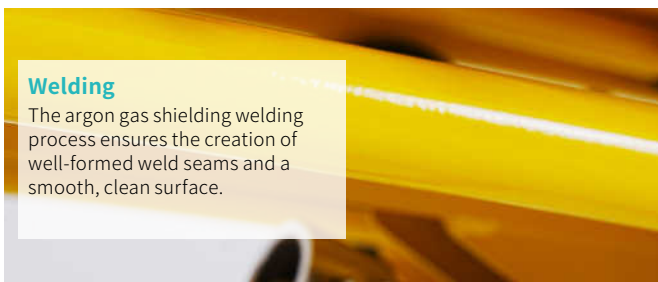


Compliant suspension system

By utilizing a compliant suspension system, different pressure points on the bed surface will not lead to the overturning of the stretcher.

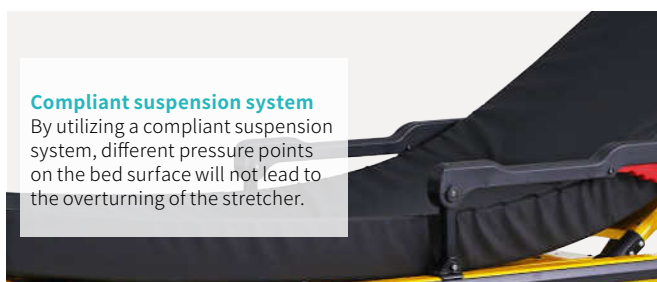

Self-adjusting system

When boarding or alighting, the self-adjusting system efficiently distributes weight, effectively preventing the common forward movement associated with stretcher loading.

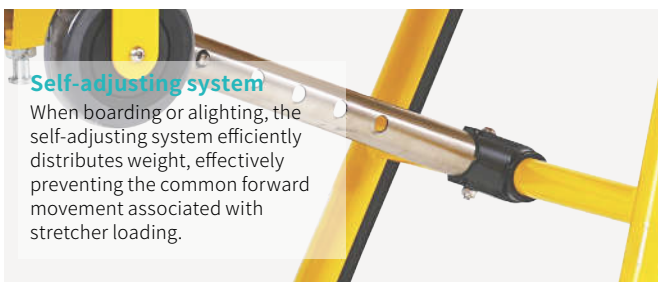

Stability

Regardless of how the patient's pressure points are distributed on the stretcher, the unique geometric design ensures absolute stability during transportation.



Ambulance platform

The stretcher's structure can be easily adjusted to match the height of any ambulance platform. Height adjustments can be effortlessly made by simply adjusting the front axle.



Folding leg mechanism

The mechanical folding leg mechanism, controlled by a rear handle, automatically folds when boarding, enabling single-person operation and easy boarding.




Low-profile chassis





A low-profile chassis structure that can move on the ground and be lifted for walking.



Side rail

Made from high-density polymer plastic, the flip-up guardrails protect the patient during transfer.

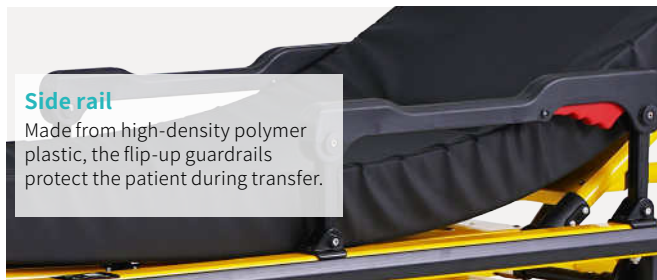

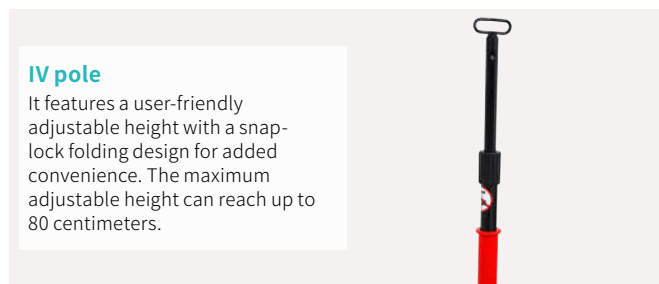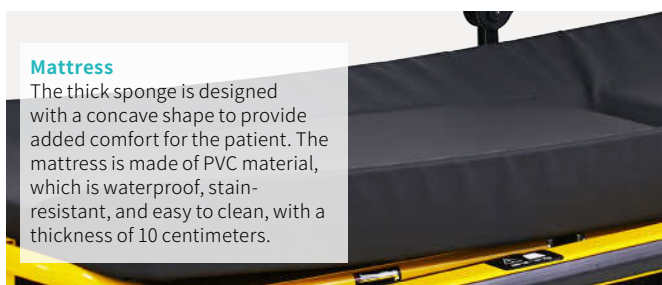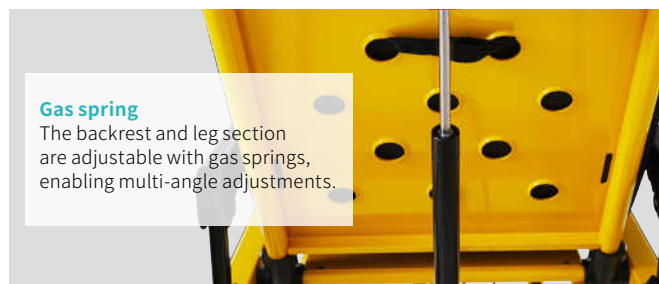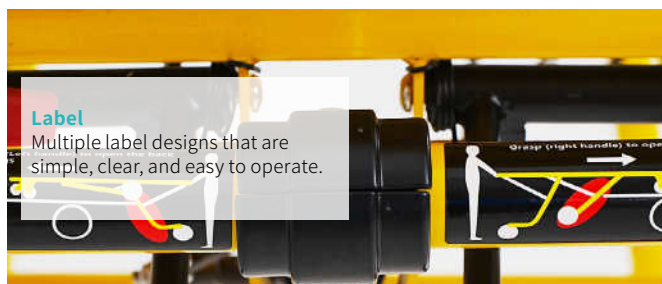

Adjustable backrest

The backboard is formed through a single molding process, allowing the head and upper body positions to be adjusted from 0° to 65°, and the foot section to be adjusted from 0° to 15°. This facilitates cleaning and disinfection and enables the adjustment of the patient's body position according to different medical requirements. It is mounted on the metal mainframe of the stretcher, below the patient mattress, to prevent secondary injuries to fracture patients during transport. When necessary, cardiopulmonary resuscitation (CPR) can be performed directly on the stretcher.


Three-stage height adjustment

The stretcher trolley features a unique three-stage height adjustment design, where each level can be easily achieved by pulling the telescopic frame. This design assists paramedics in smoothly lifting patients onto the stretcher.


Detachable stretcher

The stretcher surface is detachable and comes equipped with a locking device, ensuring secure locking during transportation onto the ambulance.



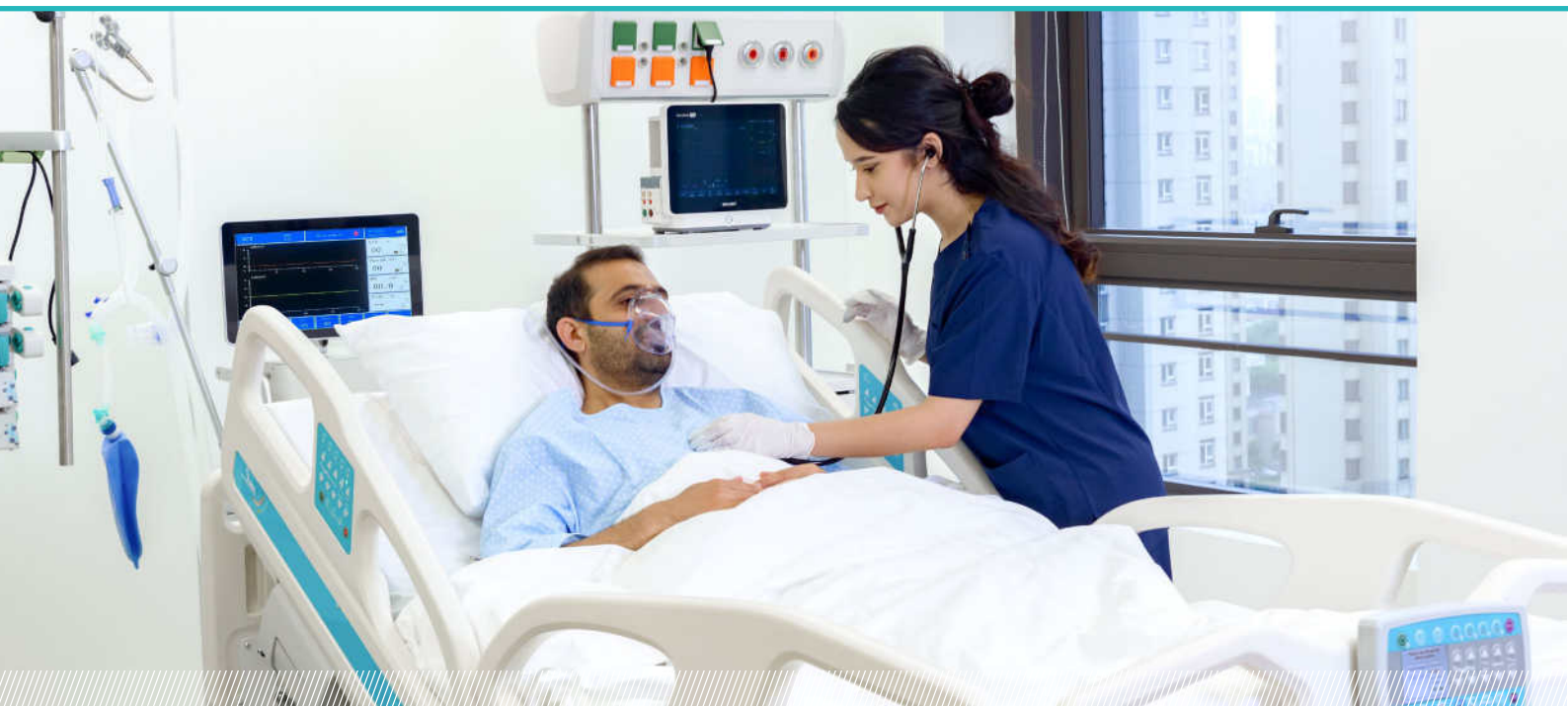

CONSTANT INNOVATION

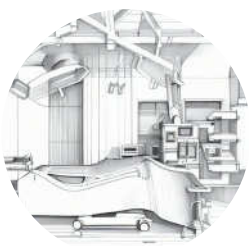


Design



Research



Manufacturing

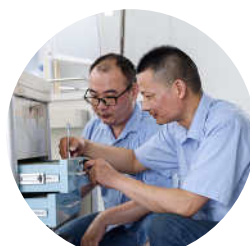

Maintenance



Service



WORLDWIDE BRANCHES

Portugal Showroom

Av. da Liberdade, nº 101
8150-101
São Brás de Alportel
Portugal
TEL: +351 289 090 212

Brazil Showroom

Rua Lauro Muller 116, Sala
1503 - Botafogo, Rio de
Janeiro - RJ, 22290-160
E-mail: mehecobrasil@
gmail.com

Indian Showroom

1460 M.I.E, Bahadurgarh,
Haryana-124507, INDIA
TEL: +91 96677 70747
E-mail: ind@saikangmedical.
com

China Headquarter

No.35 Lehong Road,
Zhangjiagang City,
Jiangsu Province, China.
TEL: +86 18021231901
E-mail: export@
saikangmedical.com

South East Asia Branch

**SAIKANG Medical
(Thailand) Co., Ltd.**
Ample Tower 909 Room 11/3
Debaratana Road,
Bangna Nuea, Bangkok
10260 Thailand
TEL: +66-936-569469

US Branch

**Healthward Technology,
INC**
17800 castleton st ste 180,
Los Angeles
TEL: +1 (380) 254 4913

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue