

CELEBRATING
ANNIVERSARY

23rd

· 2002-2025 ·



Healthier your life

C8c


ELECTRIC ICU BED

Registration code: SKD-CV

www.saikangmedical.com

SAIKANG Medical®

THE EVOLUTION OF MEDICAL FURNITURE

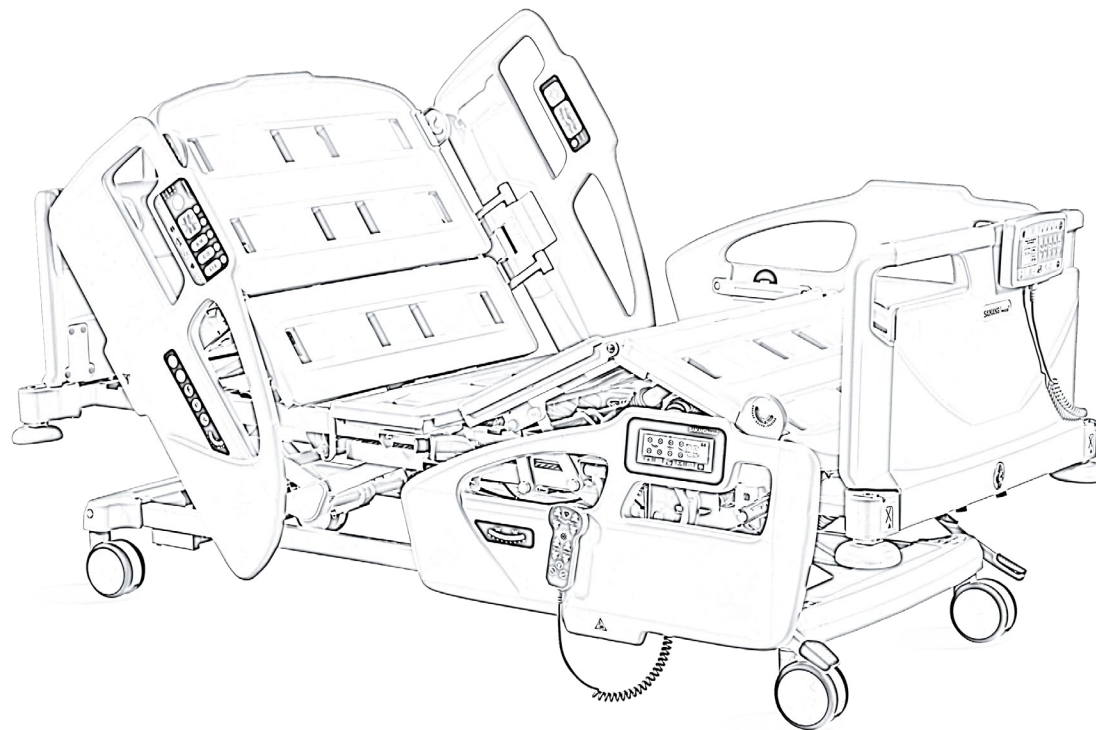

Unique aseptic design to meet hospital requirements for infection control. Seamless headboard, footboard and bed platform for easy disinfection. The width of the whole bed is 990mm, which is easy to move in the ward.

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

Panasonic Robotic ensure 360° full smooth welding.

Solid bed base made of 3.0mm thick 30×50 rectangular tube, safe working load 250kg, with a fully enclosed design for dirt protection and anti-entangling.

RACING AGAINST TIME HEALTHIER YOUR LIFE

ENGINEERED FOR PATIENT WELLNESS, DESIGNED FOR CLINICAL EXCELLENCE

Working alongside individuals in critical situations, we recognized the pivotal role of the bed.

- The question arose: to what extent can care be enhanced by improving the operator's capabilities through the bed?

- Transforming the approach, we see designing a bed tailored to the caregiver as not a constraint but an influential and adaptable instrument.


“DESIGNED TO ASSIST ME IN
WORKING MORE EFFICIENTLY”






Backrest 30-degree smart stop: The clinically recommended position for preventing respiratory complications. Research shows that the 30-degree position significantly reduces the risk of lung infections.

Auto-Regression

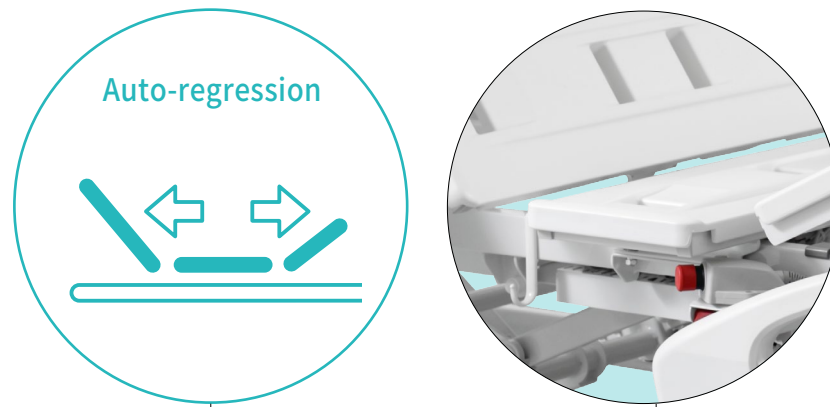

DUAL RETRACTION PATENTED DESIGN INNOVATIVE 67+57 mm

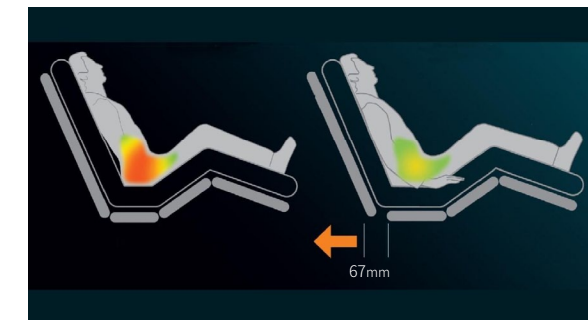

Pressure ulcers, or bedsores, result from prolonged pressure and friction on skin and underlying tissues.

- The dual retraction feature in a medical bed extends and distributes pressure evenly, reducing pressure concentration on the pelvic area.
- This design enhances patient comfort and prevents pressure ulcer occurrence proactively.
- The patented feature is applied to the backrest and thigh rest, retracting them by 67mm and 57mm respectively.
- Retracting these components extends the bed surface in the pelvic area, alleviating pressure and shear forces, thus minimizing pressure ulcer risk.
- The innovative design prioritizes patient well-being and offers a proactive solution for pressure ulcer prevention.

INTEGRATED PP BED PLATFORM

Integrated PP bed board, using advanced manufacturing technology.


5 pcs bed plates, independent installation, easy maintenance.




The back plate and the foot plate with integrated mattress retainers of 30mm height.



The detachable mattress retainers on headboard and footboard are 56mm high, 90mm wide, along with the mechanical angle indicator.

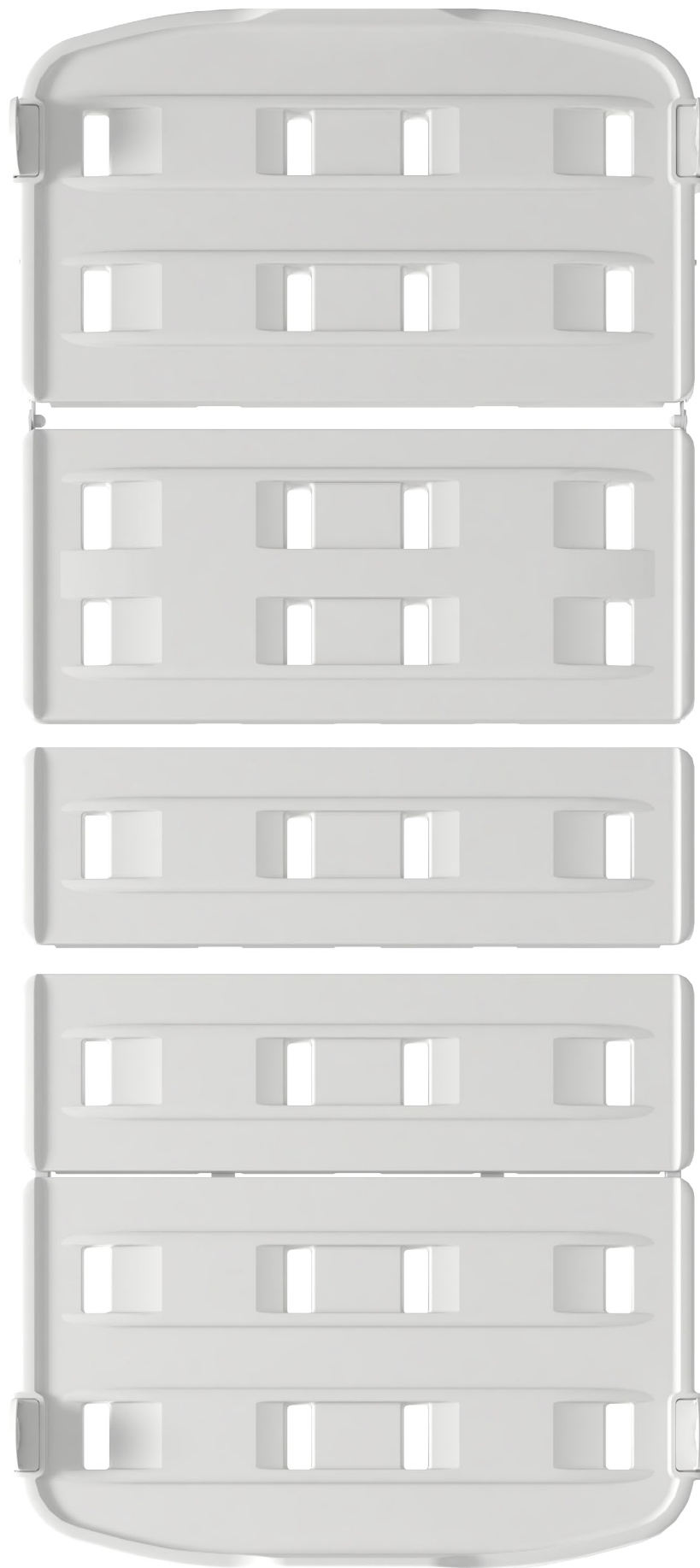


SAIKANG Medical®

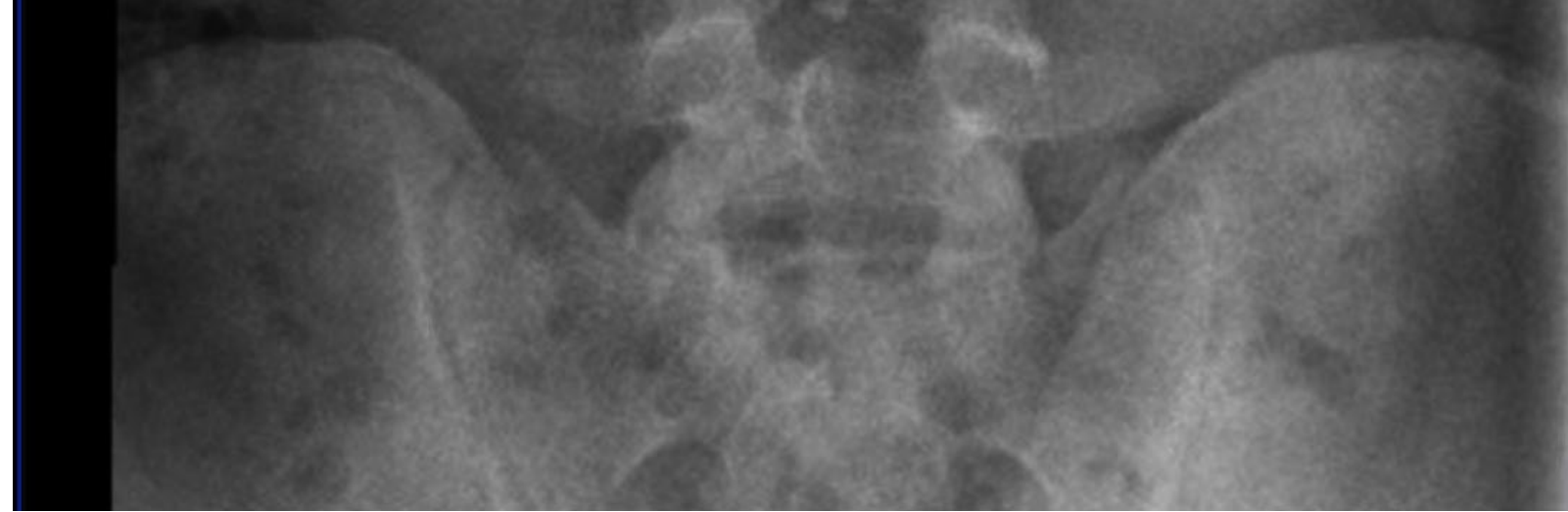

X-RAY TRANSLUCENT MATTRESS PLATFORM (BACK SECTION)

X-ray translucent mattress platform with X-ray cassette 540*365mm (Back Section).

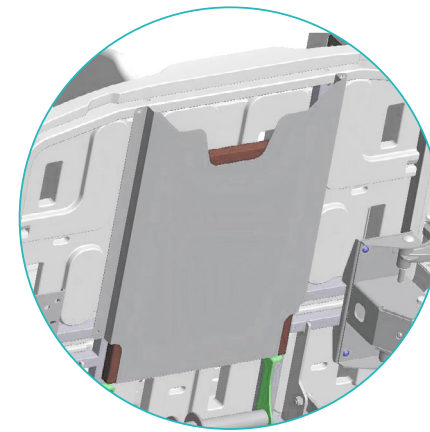



Back-rest



Bed height min.



Reverse Trendelenburg



The examination position height adjustment facilities medical staff in the examination, medication and monitoring operations, by standing up straight.



Knee-rest



Bed height max.



Trendelenburg



The CPR adjustment allows to bring the mattress platform immediately down to the horizontal.



The chair position helps prevent respiratory complications, it induces to a more physiologic position and has a psychologically positive influence.

BLOW-MOLDED SEAMLESS PP SIDERAIL



Integrated multi-function control panel, touch-screen designed for easier touch. Each touch point has light feedback to avoid misoperation; the acrylic surface is more flat, waterproof and resistant to rubbing test, convenient for patients to adjust the basic positions of the bed through the controller.



The control pannel is integrated in the siderail with the outer part having a 1.3 inch circular LCD screen to display the angle of the backrest.



Optional 6.86-inch long LCD screen on foot-end siderail to display the patient's weight.



Tuck away siderail with buffer integrated and buit-in angle indicator.



Integrated warning lights for more obvious warning.



Optional wireless charger (15W) integrated into the storage tanks on both sides of foot-end siderails.



Each siderail with a groove to fix manual controller.

HEAD & FOOT BOARD

Blow-molded seamless PP headboard and footboard.



THE BED END IS EQUIPPED WITH PROJECTION LIGHTS

The bed end is equipped with projection lights, which can respectively project the warnings;

Yellow pattern warnings: nurse call, siderail unlocked; red pattern warnings: brake unlocked, out-of-bed position; green pattern warnings: the lowest position.

These patterns are convenient for the nursing staff to observe the state of the bed and the needs of patients from outside.



NURSE CALL



SIDERAIL UNLOCKED



Plug in type, easy to install, one step unlock, knob-type unlocking mechanism for easier operation.



Crystal adhesive logo: beautiful and easy to clean.



With a limit convex point to avoid the slide of the nurse controller and hand controller.



BRAKE UNLOCKED



OUT-OF-BED POSITION



THE LOWEST POSITION

EFFORTLESS MOBILITY: SMOOTH AND SECURE REPOSITIONING WITH MINIMAL CAREGIVER STRAIN


Its compact dimensions make it ideal for use in multi-bed rooms. Good maneuverability is crucial to facilitate the work of staff and prevent collisions. The modern castors make the bed easy to move, while an optional fifth castor facilitates maneuvering around tight corners.



5-inch central locking double-sided castors. TPR tire no worn out after running 10 km.

Anti-winding hard shell cover, united forming without bolts.

Pass dynamic test: Successfully pass 40,000 obstacle tests without deformation, under the standard 200kg load of the entire bed.

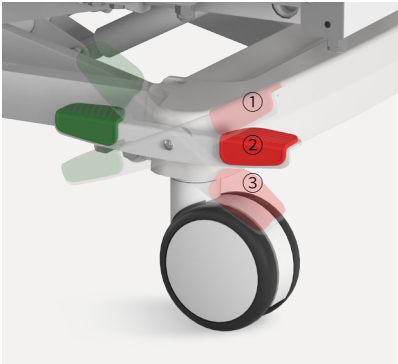

- ①-speed steering,
- ②-speed free,
- ③-speed locking.



PATIENT WEIGHT UNDER CONTROL AT ALL TIMES



The exclusive weighing system was designed to guarantee a high level of precision independently from the mattress platform's position.

Thanks to the 4 sensors located at the corners of the bed, the load is better distributed and the reading level more accurate.

The weighing data can be recorded manually whenever necessary or automatically every 24 hours.

While the weight of the patient is monitored it is possible to set up different types of alarm: exceeded threshold weight alarm or "bed exit" alarm.



The bed platform can be extended by 20 cm, providing added comfort for taller patients.



Back-up Battery
Battery – Can operate 18 cycles of lifting and lowering with 100kg load bearing, and up to 130 cycles without load bearing.

QUICK EMERGENCY ADJUSTMENTS

During resuscitation, it serves as a crucial aid. The shock position can be instantly activated using the easily accessible emergency lever at the head end or by using the control devices. Quickly removable headboards and footboards allow unimpeded access to the patient.



3.5 inch LCD screen, flexible spiral cable for easy mobility and operation in optimal position. The nurse can control all necessary functions.

One 5-function hand controller (backrest, legrest, combined back and leg adjustment, bed height adjustment, Trendelenburg, reverse Trendelenburg, CPR, button backlight, flashlight, power switch) – self-developed system.



The CPR position can be set even faster using the mechanical emergency lever under the mattress base.



One button stop, prevent the maloperation.



LINAK motor, with UL/ROHS/EN standard. Noise under 50dB, with lifetime 20,000 times movement.

Night Light

It is convenient for operators and caregivers to observe the situation under the bed even when the light is insufficient and it is energy-saving and environmentally friendly, and the induction switch does't need manual operation.

Product Information




IV Pole

The IV pole is made of stainless steel and has four plastic hooks that can bear up to 15 kg.



Drainage Hooks

Plastic drainage hook for flexible position changing.

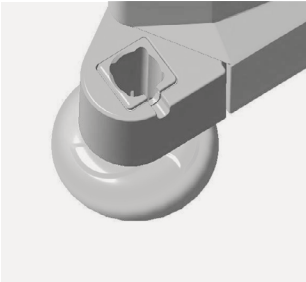

Trapeze Bracket

Integrated design for both IV pole and Lifting pole.



Mechanical Lever

Footboard with a mechanical lever for calf raise, offering 10 height levels to achieve a semi-Fowler position to improve the breathing of shock patients.

Back-rest adjustment

Knee-rest adjustment

Height adjustment

Reverse trendelenburg

Trendelenburg

Auto-regression

Position line

Angle indicator

Instant release

Cardiac Chair Position

Built-in LED torch

Weighing system

Bed extension

LINAK Motor

Nurse Call

Siderail Unlocked

Brake Unlocked

Lowest Position


Out-of-bed Position Warnings

Model	C8c
External size (LxW)	2180x990mm±10mm
Mattress platform	1995x860mm±10mm
Safe working load	250kg
Height adjustment	400-850mm±10mm
Back-rest adjustment	0-70° (±10°)
Knee-rest adjustment	0-28° (±10°)
Reverse trendelenburg	0-17° (±2°)
Trendelenburg	0-17° (±2°)
Angle indicator	0-90°
Auto-regression	67+57mm
Position line	30°
Bed Extension	200mm
LINAK Motor	●
5" Double Side Castors with directional wheel function	●
PP Removable C-series Headboard and Footboard Φ113 Bumper Wheels	●
PP 4-pcs C-Series Siderails	●
Head-end Siderails with Integrated Controller and Warning Lights	●
Foot-end Siderails with Glasses Holder, LCD Screen and 15W Wireless Charger	●
PP Detachable Bed Platform With Integrated Mattress Retainer	●
Mechanical Lever on Leg Board for Calf Raise	●
Backrest with X-ray Cassette Holder	●
Patent Angle Indicator for Back-rest and Trendelenburg Position	●
IV Pole 4 Hooks	●
IV Pole Prevision	●
Drainage Hook	●
Φ32mm Monky Pole Holes	●
Hand Cotroller	●
Nurse Cotroller	●
Cardiac Chair Position	●
Electric CPR	●
Manual CPR	●
Auto Break at 30°	●
Emergency Stop Button	●
Night Light	●
Bed Extension	●
Weighing System	●
Bed End with Projection Lights	●
Nurse Call	●
Siderail Unlocked	●
Brake Unlocked	●
Lowest Position	●
Out-of-bed Position Warnings	●

Note: ●-YES, ○-Optional



CONSTANT INNOVATION

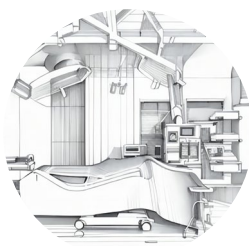


Design



Research



Manufacturing

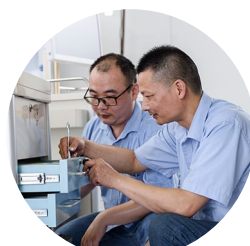

Maintenance



Service

WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited
Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC
I7800 CASTLETON ST STE 180.
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
 +86 18021231901 ·  24/7
 export@saikangmedical.com
 www.saikangmedical.com



Company video



General catalogue