

CELEBRATING
ANNIVERSARY

23th

· 2002-2025 ·



Healthier your life

X9x


ELECTRIC ICU BED

Registration code: SKD-CV

www.saikangmedical.com

SAIKANG Medical®

THE EVOLUTION OF
MEDICAL FURNITURE

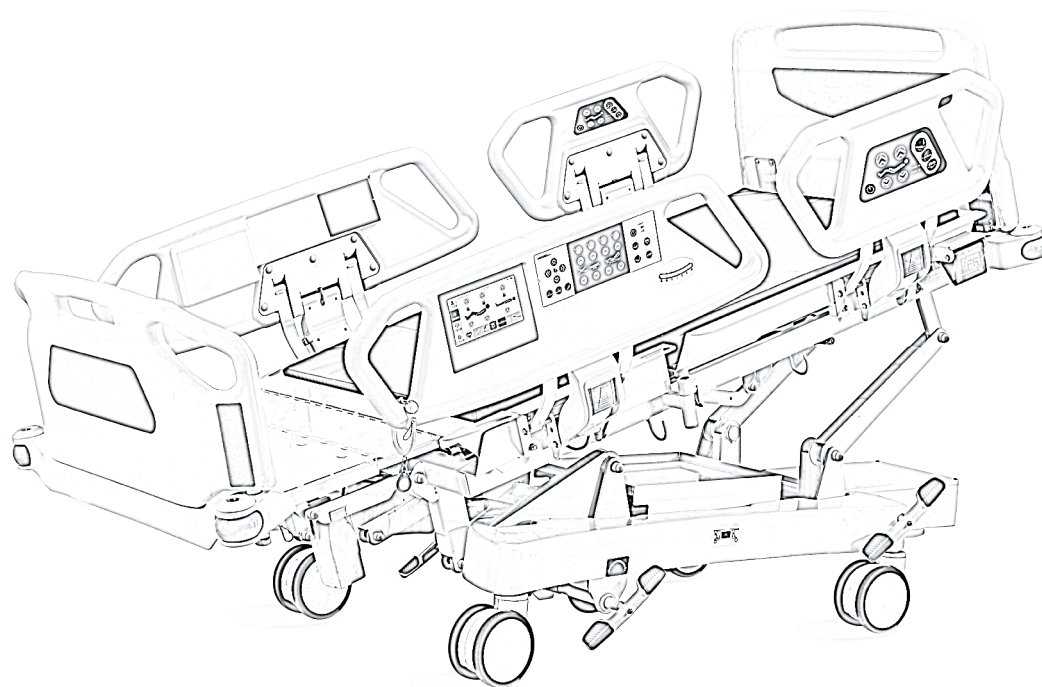

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

Panasonic Robotic ensure 360° full smooth welding.

RACING AGAINST TIME
HEALTHIER YOUR LIFE

ENGINEERED FOR
PATIENT WELLNESS,
DESIGNED FOR CLINICAL
EXCELLENCE

Working alongside individuals in critical situations, we recognized the pivotal role of the bed.

- The question arose: to what extent can care be enhanced by improving the operator's capabilities through the bed?

- Transforming the approach, we see designing a bed tailored to the caregiver as not a constraint but an influential and adaptable instrument.


“DESIGNED TO ASSIST ME IN
WORKING MORE EFFICIENTLY”






Linak motor from Denmark,
ensure durable stability. With
UL/ROHS/EN Standard.

Backrest Regression

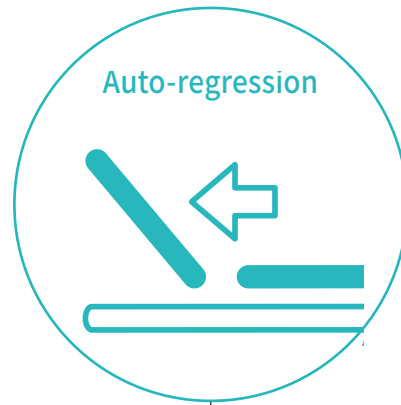

BACKREST RETRACTION PATENTED DESIGN INNOVATIVE 75mm

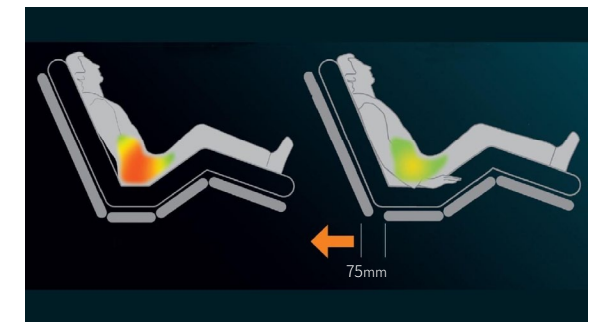

Pressure ulcers, or bedsores, result from prolonged pressure and friction on skin and underlying tissues.

- Auto-Regression 75 mm $\pm$ 10mm

- When repositioning the backrest , the back plate retreat 75mm to minimises the pressure and shear forces in the patient's abdomen region. This results in a significant reduction in the tugging on the patient's wound.

- This design enhances patient comfort and prevents pressure ulcer occurrence proactively.

- Retracting these components extends the bed surface in the pelvic area, alleviating pressure and shear forces, thus minimizing pressure ulcer risk.

- The innovative design prioritizes patient well-being and offers a proactive solution for pressure ulcer prevention.

SOLID YET LIGHT MATTRESS PLATFORM

8 Functions:

Backrest adjustment,
Kneerest adjustment,
Height adjustment,
Trendelenburg,
Reversed trendelenburg,
Bed extension,
Cardiac chair position,
Zero position.



7cm PU mattress



X-RAY EXAMINATION WITHOUT MOVING THE PATIENT

Wide area of radiotransparency. The Phenolic Resin Back Platform has an X-Ray cassette holder for chest x-ray examinations, without moving the patient and at any angle of the backrest.

X-ray translucent mattress platform (Back section) with X-ray cassette (size 550*378mm)

The X-ray tray holder fixed at the back section, easy for X-ray cassette film horizontally put out & put in, to avoid patient movement.





Back-rest



Bed height min.



Reverse Trendelenburg



Knee-rest



Bed height max.



Trendelenburg




One button nursing operating position: bed platform 780mm high from the ground, helping the nurse to give fast treatment in a very comfortable height.



Electric CPR button is setted inside control panel and nurse control panel, one button to make the back plate flat.



The position of the cardiac chair can be easily adjusted with a single button, optimizing respiratory function and cardiac output.

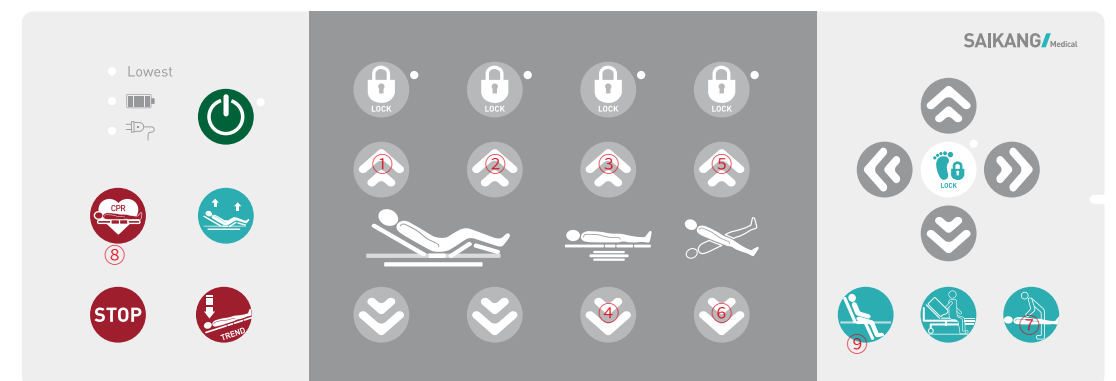

SAFE AND EASY TO OPERATE

ADVANTAGES FOR NURSING CARE

Just one-hand operation required to effortlessly release and lower the split safety siderail.

Tuck away side rail with buffer.

Built-in controller and angle indicator.



7-inch touch screen, instead of the original button screen, more convenient for operation

Screen can store the patient and doctor's information, users manuals, weight record.



CUSTOMISED CONTROL CONCEPTS

Nursing staff and patients require safe and convenient bed operation options. Our 'light' version control panel provides customized features both inside and outside the bed. By prioritizing essential functions, operation becomes significantly simpler. An easily accessible handset with ergonomic buttons is also provided as an alternative.

Easy operation using the control panel “light”

Staff and patients benefit from the tailor-made convenience and care functions that are available on either side of the "light"-version control panel. The panel can be integrated in the split safety sides if required.



Easy bed adjustment with the “control panel light” version

The “control panel light” version is always within easy reach of patients. The straightforward adjustment options are easy to understand even for people with cognitive impairment.

The “light”-version control panel is always at eye level for the patient.


ALARM EXIT SYSTEM



Used to monitor the movement of patient during transitional rehabilitation.



The lowest height position minimises the rate of fall in the event of any fall. This significantly reduces the risk of injury.

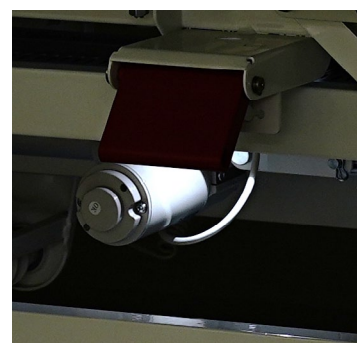

Night light underneath the bed helps patient getting up at night.

QUICK EMERGENCY ADJUSTMENTS

During resuscitation, it serves as a crucial aid. The shock position can be instantly activated using the easily accessible emergency lever at the head end or by using the control devices. Quickly removable headboards and footboards allow unimpeded access to the patient.

Nurse lock function, one button flat.



Electric CPR—one key to release of backpart

Dural side handle to release of backpart with damping device which can let the back rest fall down slowly to ensure the safety of the patient, single hand operation.




The CPR position can be set even faster using the mechanical emergency lever under the mattress base.



The headboard and footboard can be unlocked and removed without tools to allow improved access to the patient.

PP bed ends One Step Removable, meet the first-aid requirements of high-end wards and ICU room

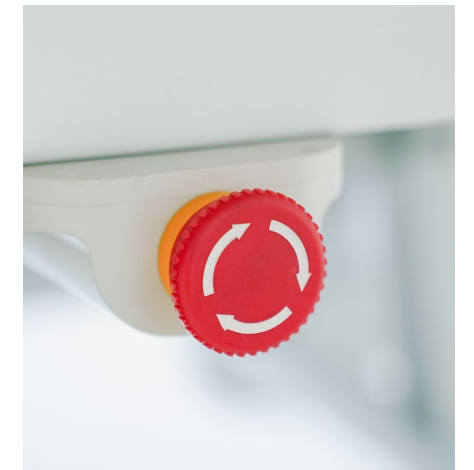

One-button emergency stop switch can quickly stop all going movement


EFFORTLESS MOBILITY: SMOOTH AND SECURE REPOSITIONING WITH MINIMAL CAREGIVER STRAIN

Its compact dimensions make it ideal for use in multi-bed rooms. Good maneuverability is crucial to facilitate the work of staff and prevent collisions. The modern castors make the bed easy to move, while fifth castor facilitates maneuvering around tight corners.



6" diameter central locking castors with TPR tire, no worn-out after running 30KM, united forming without bolts

Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time

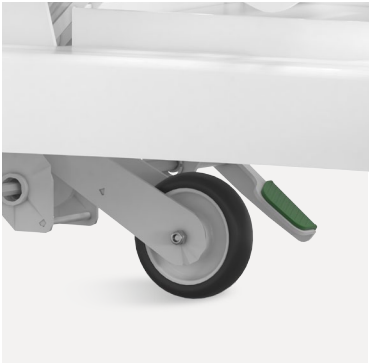

The 5th directional wheels, during the movement can well control the direction.

PATIENT WEIGHT UNDER CONTROL AT ALL TIMES




The exclusive weighing system was designed to guarantee a high level of precision independently from the mattress platform's position.


Thanks to the 4 sensors located at the corners of the bed, the load is better distributed and the reading level more accurate.

The weighing data can be recorded manually whenever necessary or automatically every 24 hours.

Weighing system, accurate to 0.01Kg, with data storage system store 10times datas.

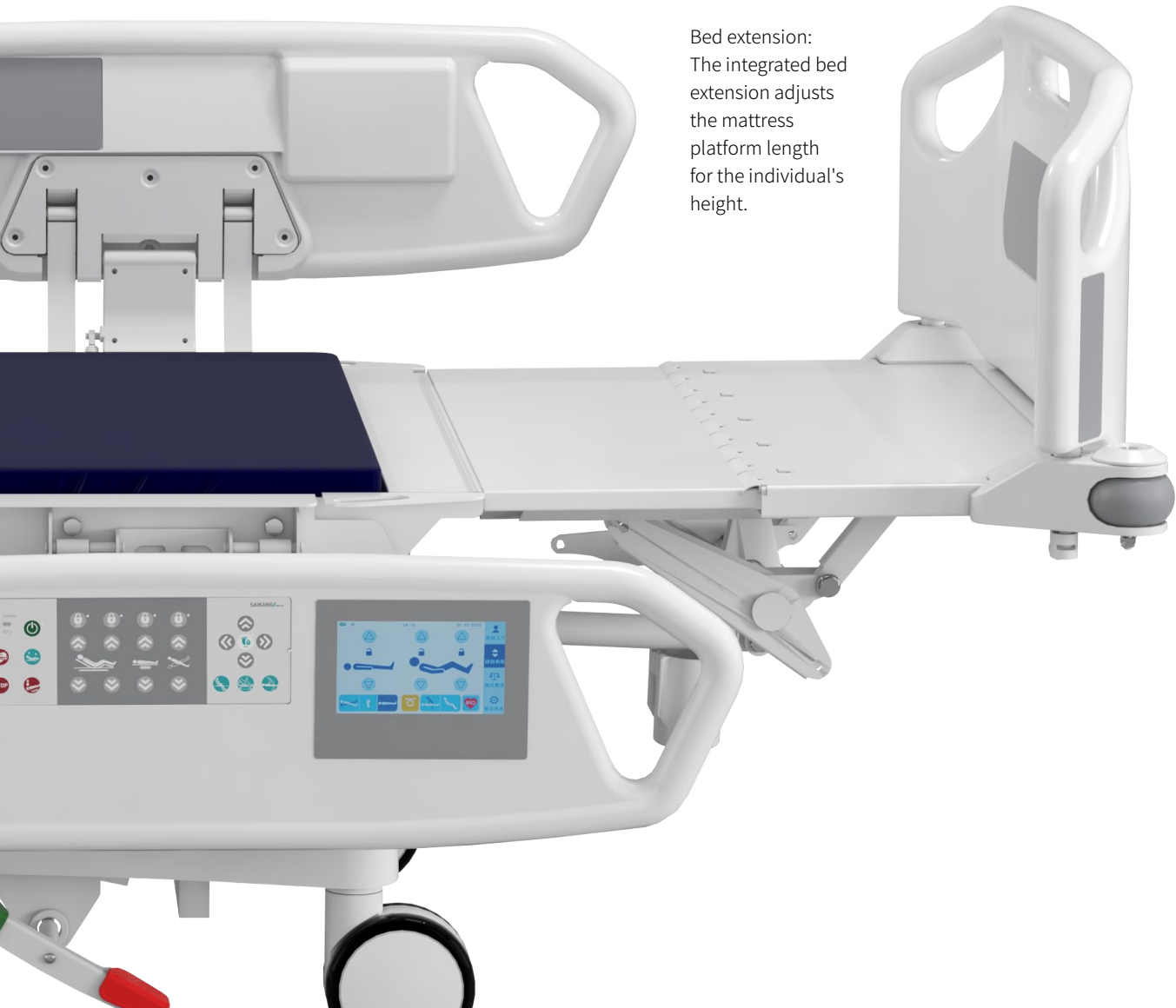


Bed extension:
The integrated bed extension adjusts the mattress platform length for the individual's height.

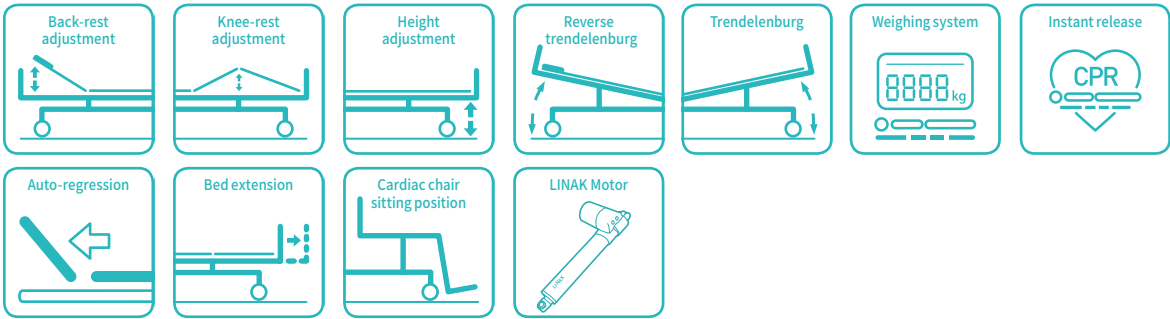

The I.V. pole is made of stainless steel and has four plastic hooks that can bear up to 15 kg of weight.



Drainage hook 4pcs.

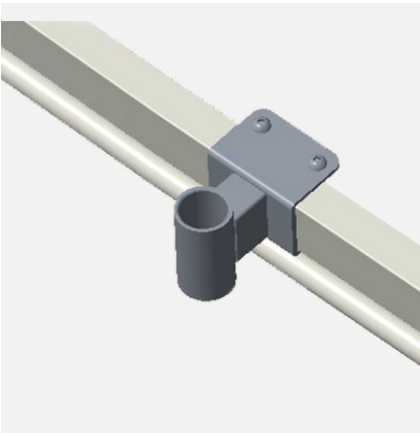

Trapeze Bracket.

Options



Trapeze Pole
32mm dia heavy duty trapeze pole.



Back up battery, 18 times movement with 100kg bearing, 130 times movement without bearing.


Product Information

Model	X9x
External size (LxW)	(2120-2320)×1010 mm
Mattress platform	(1750+200)×830 mm
Safe working load	250 kg
Height adjustment	550-950 mm
Back-rest adjustment	0-75° (±10°)
Knee-rest adjustment	0-20° (±10°)
Reverse trendelenburg	0-18° (±2°)
Trendelenburg	0-18° (±2°)
Leg extension board Adjustment	0-68° (±2°)
Angle indicator	0-15°
Backrest-regression	75 mm
LINAK Motor	●
6" Double Side Secure®Castors, Central Locking	●
ABS Side Rail with Integrated Control Panel	●
ABS Bed Ends One Step Removable	●
Φ85mm Wall Bumper Wheel	●
IV Pole	●
IV Pole Prevision	●
Drainage Hook	●
5th Directional Wheel Pedal	●
Central Locking Pedal	●
X-ray Translucent Mattress Platform(Backrest)	●
Back-up Battery	●
Night Light	●
Trapeze Socket	●
7cm Mattress	●
Bed Push Handle for Emergency Transfer	●
Emergency Stop Button	●
Weighing system	●
Guardrail LCD Screen with sensor key	●
Manual CPR	●
Electric CPR	●
Monkey Pole	○

Note: ●-YES, ○-Optional



CONSTANT INNOVATION

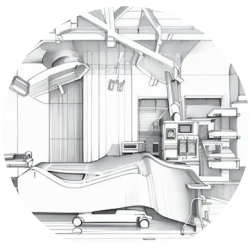


Design



Research



Manufacturing

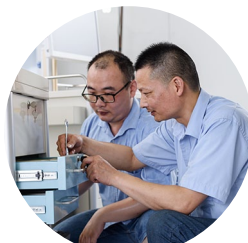

Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited
Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC
17800 CASTLETON ST STE 180,
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue