

CELEBRATING
ANNIVERSARY

23rd

· 2002-2025 ·



Healthier your life

Y8y8t


ELECTRIC ICU BED

Registration code: SKD-CV

www.saikangmedical.com

SAIKANG Medical®

THE EVOLUTION OF MEDICAL FURNITURE

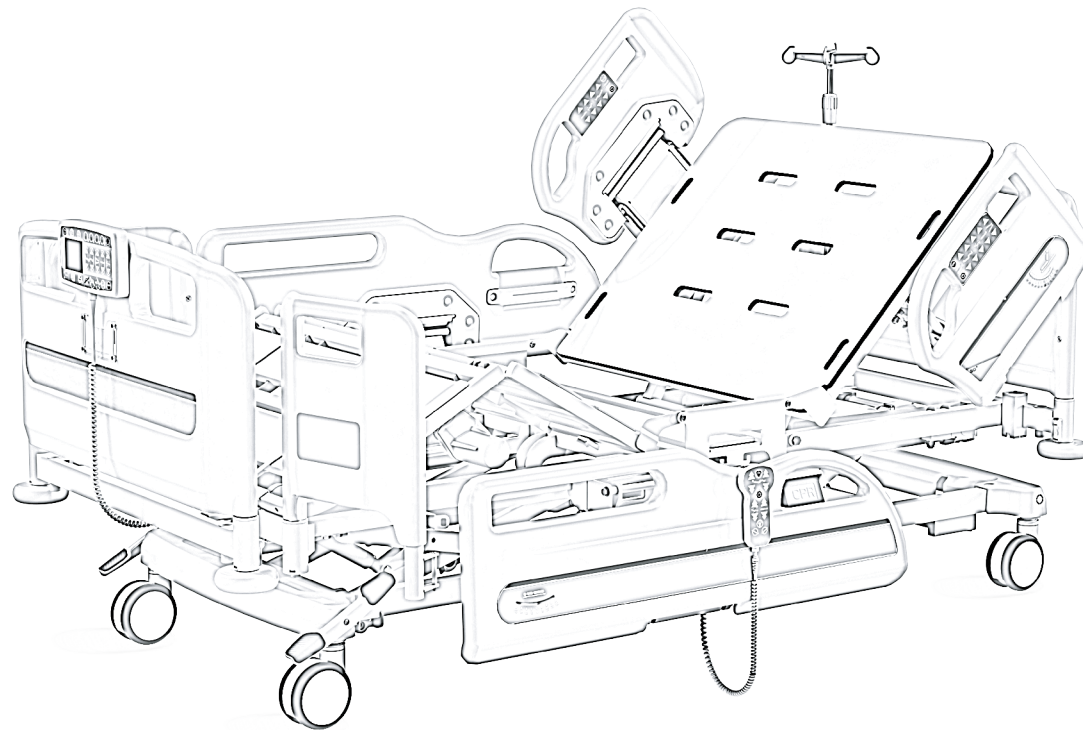

The unique germ-free design meets hospitals' strict infection control requirement. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization.

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

Panasonic Robotic ensure 360° full smooth welding.

RACING AGAINST TIME HEALTHIER YOUR LIFE

ENGINEERED FOR PATIENT WELLNESS, DESIGNED FOR CLINICAL EXCELLENCE


Working alongside individuals in critical situations, we recognized the pivotal role of the bed.

- The question arose: to what extent can care be enhanced by improving the operator's capabilities through the bed?

- Transforming the approach, we see designing a bed tailored to the caregiver as not a constraint but an influential and adaptable instrument.

“DESIGNED TO ASSIST ME IN
WORKING MORE EFFICIENTLY”





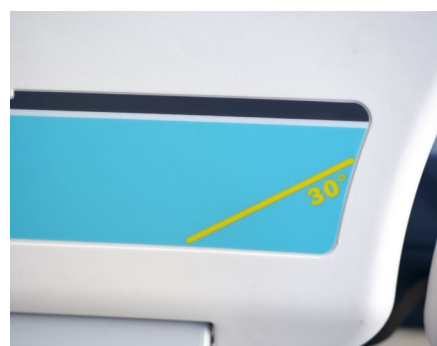


Backrest 30-degree smart stop: The clinically recommended position for preventing respiratory complications. Research shows that the 30-degree position significantly reduces the risk of lung infections.

Auto-Regression

Double auto-regression

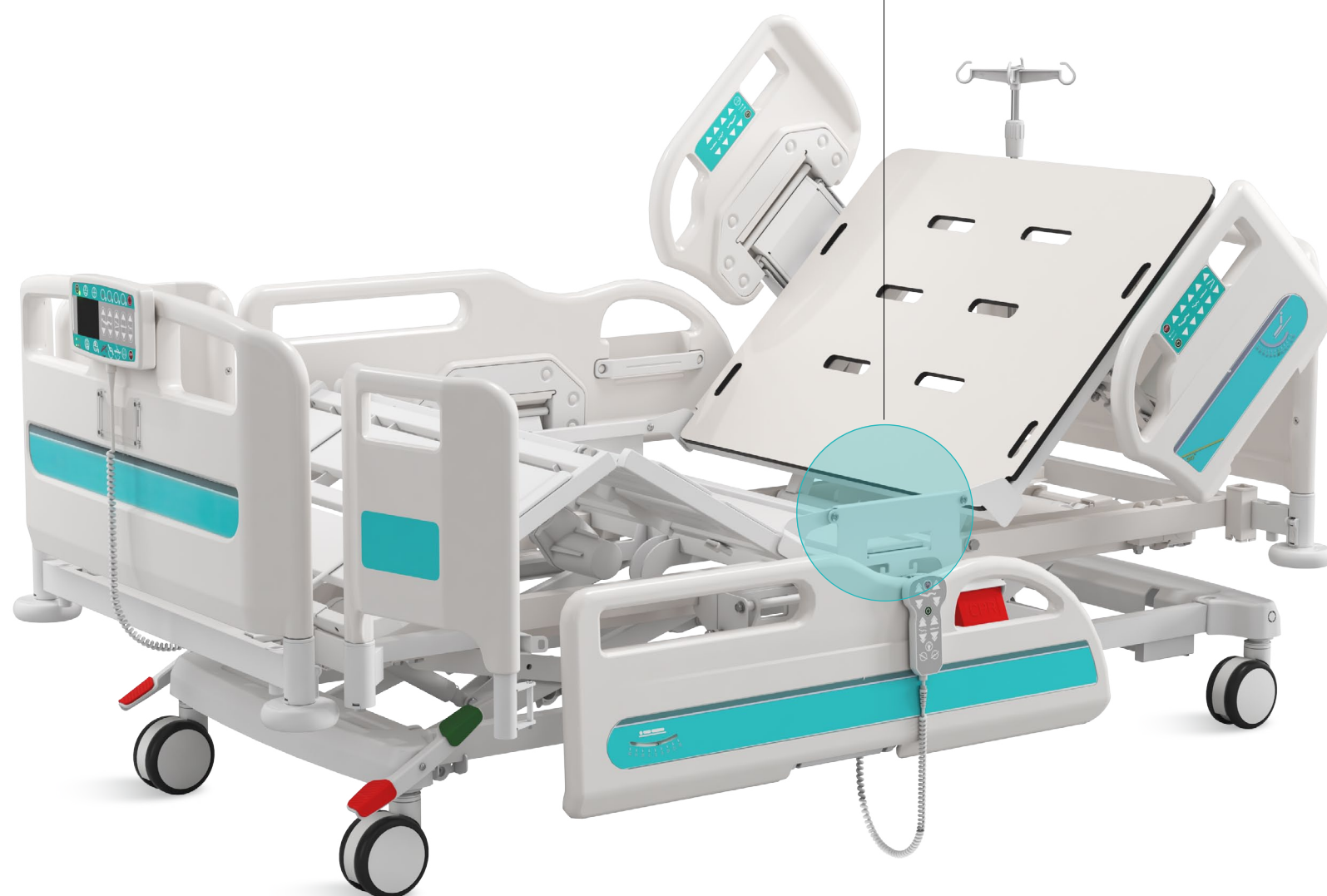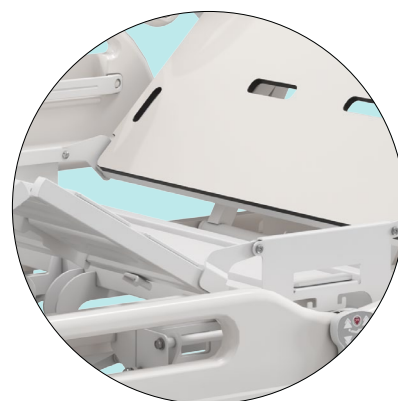

DUAL RETRACTION PATENTED DESIGN INNOVATIVE 70+40 mm


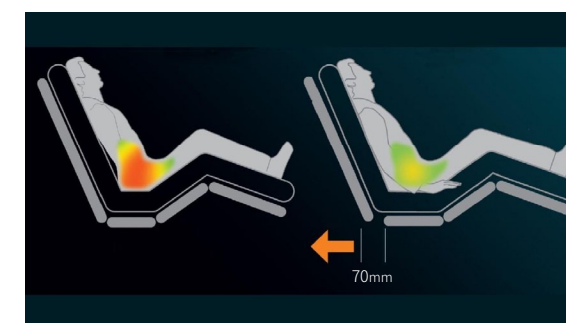


Pressure ulcers, or bedsores, result from prolonged pressure and friction on skin and underlying tissues.

- The dual retraction feature in a medical bed extends and distributes pressure evenly, reducing pressure concentration on the pelvic area.
- This design enhances patient comfort and prevents pressure ulcer occurrence proactively.
- The patented feature is applied to the backrest and thigh rest, retracting them by 70mm and 40mm respectively.
- Retracting these components extends the bed surface in the pelvic area, alleviating pressure and shear forces, thus minimizing pressure ulcer risk.
- The innovative design prioritizes patient well-being and offers a proactive solution for pressure ulcer prevention.

JOINT-FREE TECHNOLOGY

5 pieces of removable blow moulded boards to minimise the slots where unwanted infections may hide.

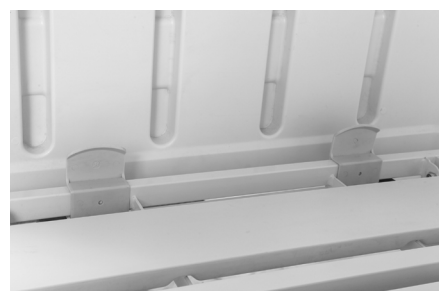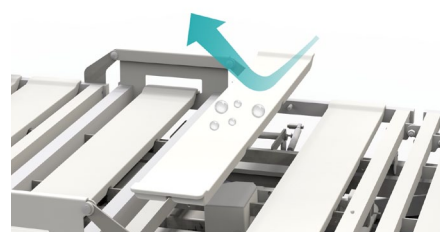

With 2 pieces of 52*70mm plastic mattress retainers.



X-RAY EXAMINATION WITHOUT MOVING THE PATIENT

Wide area of radiotransparency. The Phenolic Resin Back Platform has an X-Ray cassette holder for chest x-ray examinations, without moving the patient and at any angle of the backrest.

X-ray translucent mattress platform (Back section) with X-ray cassette (size 550*378mm)

The X-ray tray holder fixed at the back section, easy for X-ray cassette film horizontally put out & put in, to avoid patient movement.



SAIKANG Medical®



Back-rest



Bed height min.



Reverse Trendelenburg



One button nursing operating position:
bed platform 780mm high from the
ground, helping the nurse to give fast
treatment in a very comfortable height.



Knee-rest




Bed height max.



Trendelenburg



Electric CPR button is setted inside
control panel and nurse control panel,
one button to make the back plate flat.

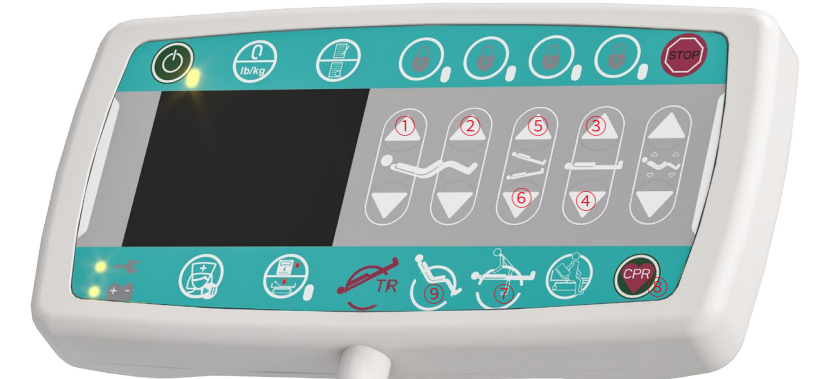

The position of the cardiac
chair can be easily adjusted
with a single button,
optimizing respiratory
function and cardiac
output.



SAFETY SIDE RAILS

The designed guardrail fully considered to provide enough space and hand grasp for patient to get out of bed, and provide extremely high safety.

Two steps unlocking system, to prevent patients from accidentally lowering the guardrail (patented design).



Glasses holder



Angle indicator



Backrest 30° indicator line



Backrest angle display at 30° and 45°




Just one-hand operation required to effortlessly release and lower the split safety siderail.

All it takes for the patient to get in and out of bed unhindered is to lower one safety siderail element in the centre of the bed.



5th and 6th extra side rail, for 100% fully protection.

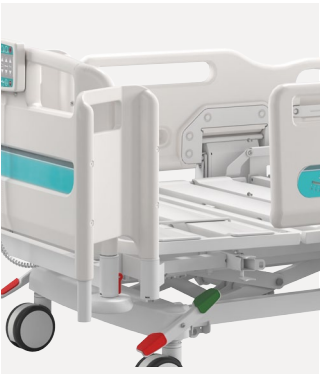

CUSTOMISED CONTROL CONCEPTS

Nursing staff and patients require safe and convenient bed operation options. Our 'light' version control panel provides customized features both inside and outside the bed. By prioritizing essential functions, operation becomes significantly simpler. An easily accessible handset with ergonomic buttons is also provided as an alternative.

Easy operation using the control panel “light”

Control panel in the siderails: Numerous studies have shown that the built-in controller located in the head guardrail is easier to reach, allowing the patient to adjust the basic position of the bed through the panel.



Easy bed adjustment with the “control panel light” version

The “control panel light” version is always within easy reach of patients. The straightforward adjustment options are easy to understand even for people with cognitive impairment.

The “light” -version control panel is always at eye level for the patient.

HAND CONTROLLER AND NUSE CONTROLLER



NURSE CONTROLS WITHIN REACH

Clear and easy to use controllers which is at the patient's disposal

QUICK EMERGENCY ADJUSTMENTS

During resuscitation, it serves as a crucial aid. The shock position can be instantly activated using the easily accessible emergency lever at the head end or by using the control devices. Quickly removable headboards and footboards allow unimpeded access to the patient.

Follow EN60601 Safety Standard



The orange emergency buttons on the outer side of the "light"-version control panel enable staff to act quickly without taking their eyes off the patient.



The CPR position can be set even faster using the mechanical emergency lever under the mattress base.



Satellite control panel for nurse: flexible spiral cable for easy mobility and operation in optimal position. The nurse can control all standard functions and automatic chair position.



Hand controller for patients and nurses, ergonomic design.




Built-in LED torch, backlit keyboard design for night use.



The headboard and footboard can be unlocked and removed without tools to allow improved access to the patient.


Patented design Detachable head & foot board, using the eccentric locking structure, one step quick removable

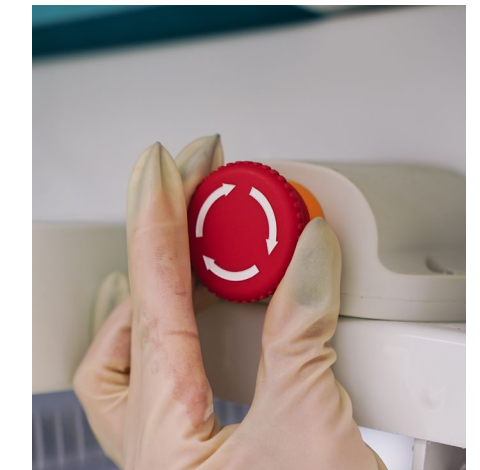

One-button emergency stop switch can quickly stop all going movement

EFFORTLESS MOBILITY: SMOOTH AND SECURE REPOSITIONING WITH MINIMAL CAREGIVER STRAIN

5" diameter central locking castors, three-position brake with directional function.
TPR tire, no worn-out after running 30KM, united forming without bolts
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time
1-speed unlock; 1-speed universal, 1-speed directional.

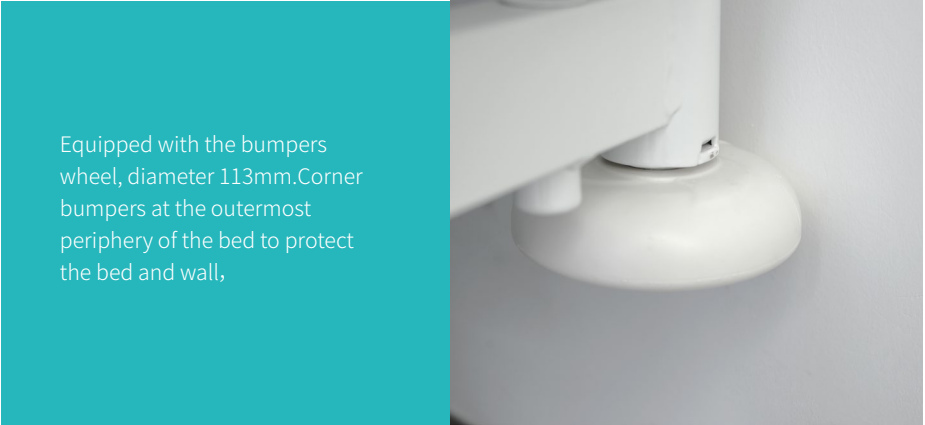

The double-folded large-area rubber design, horizontal and vertical two-sided rubber, allows the nurse to touch the brakes maximally, saving effort.




SAFE HYGIENE

Easy manual reprocessing

During a pandemic, safe and verifiable hygiene is more crucial than ever. Its low-recess design simplifies manual reprocessing for the hygiene team, saving time, physical effort, and personnel resources.



No inaccessible recesses in the headboard or footboard and safety sides..

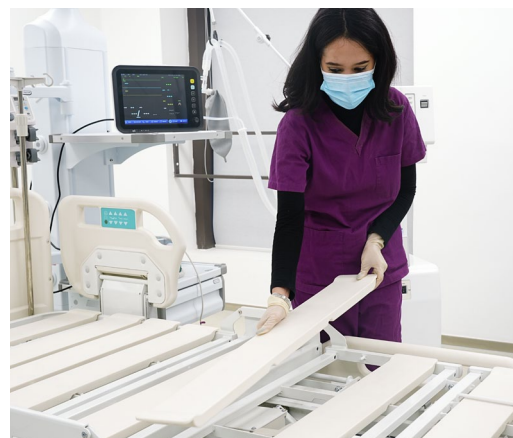

Level mattress base with easy to remove plastic cover.



Tidy cable routing.

FALL PREVENTION AT NIGHT

Patients with physical or cognitive limitations should be protected from fall injuries at night – even without the use of liberty depriving measures. The low position considerably reduces the risk in the event of a fall.




The lowest height position minimises the rate of fall in the event of any fall. This significantly reduces the risk of injury.



The under bed light helps patients by lighting their way around the room at night.



The I.V. pole is made of stainless steel and has four plastic hooks that can bear up to 15 kg of weight.

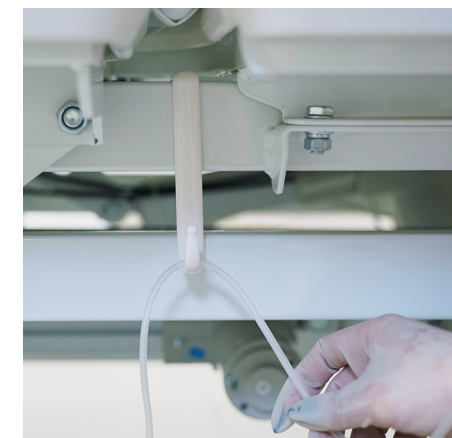

2 pieces of plastic drainage hook



4 pieces of drainage hook integrated with the bed platform.

PATIENT WEIGHT UNDER CONTROL AT ALL TIMES



The exclusive weighing system was designed to guarantee a high level of precision independently from the mattress platform's position.


Thanks to the 4 sensors located at the corners of the bed, the load is better distributed and the reading level more accurate.

The weighing data can be recorded manually whenever necessary or automatically every 24 hours.

While the weight of the patient is monitored it is possible to set up different types of alarm: exceeded threshold weight alarm or "bed exit" alarm.



The bed platform can be extended by 20 cm, providing added comfort for taller patients.

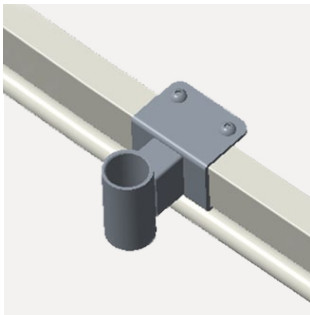

Trapeze Bracket




Back up battery, 18 times movement with 100kg bearing, 130 times movement without bearing.

SAFETY STANDARD

Follow EN60601 Safety Standard

Motor protection

The elevators are equipped with electric protection against overload, therefore the bed will automatically stop when an elevator is overloaded. Once the overload is removed, the bed can continue to work without interference.

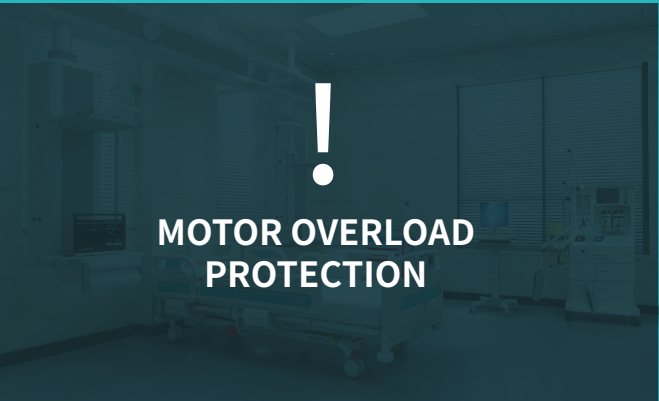

Trusted electric system

LINAK motor, with UL/ROHS/EN approval. Noise under 50dB, with lifetime 20,000times movement.




Go button

Use the GO button to activate the bed; protects the patient against an injury caused by any unintentional movement of the bed. Control system is automatically locked after 3 min.



Lowest position indicator


The lowest position indicator design allows medical staff to easily monitor the bed's safe position, even at night or from a distance. A green light indicates when the bed is at its lowest height.

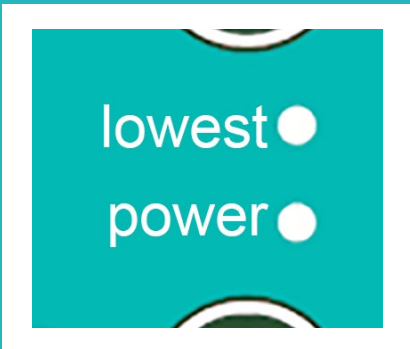


The one-touch bed exit feature

Along with the elevating backrest and the dual holes in the backrest guardrails, effectively support patients in standing up and assist them in safely getting out of bed.



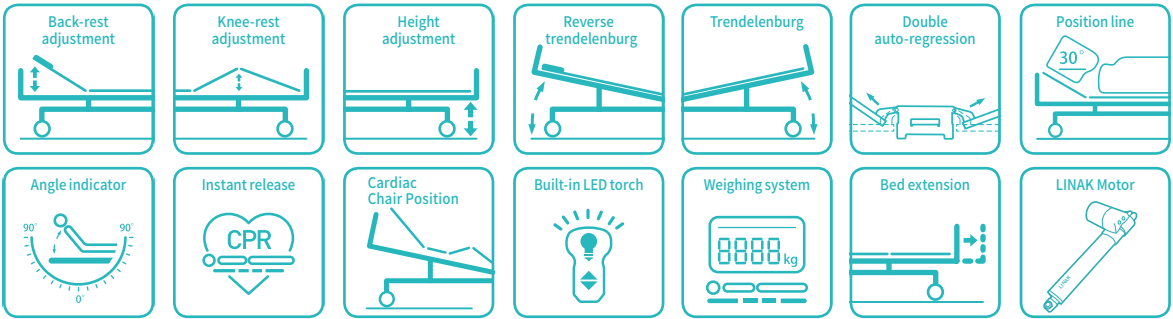

Options

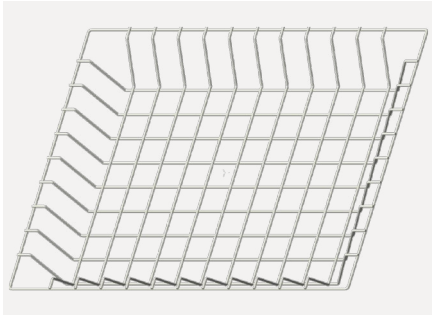


Utility Shelf

Detachable unitily shelf for patients' personal stuff.




Trapeze Pole

32mm dia heavy duty trapeze pole.

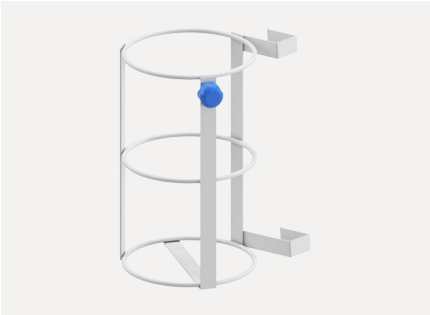

Oxygen Cylinder Holder

Optional color (MOQ100pcs- color changing)

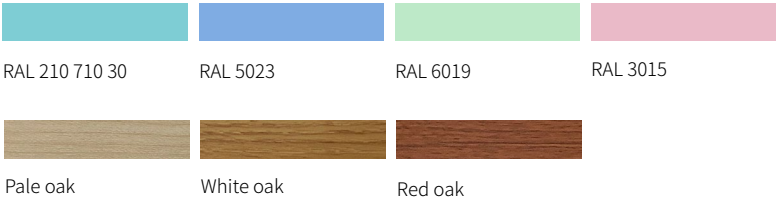


Product Information

Model	Y8y8t
External size (LxW)	2150x1060 mm
Mattress platform	1925x880 mm
Safe working load	250 kg
Height adjustment	390-840 mm
Back-rest adjustment	0-75° (±10°)
Knee-rest adjustment	0-33° (±10°)
Reverse trendelenburg	0-17° (±2°)
Trendelenburg	0-17° (±2°)
Angle indicator	0-90°
Auto-regression	70+40 mm
Position line	30°
LINAK Motor	●
5" Double Side Castors, with three-grade caster pedal	●
PP Oasis Side Rail Extra Height, with nurse control panel	●
PP Oasis Bed Ends with Safe Lock	●
PP Bed Platform,Removable&Washable	●
X-ray Transparent Phenolic Resin Back Platform	●
Plastic Socket Inside PP Guardrail for Phone or Glasses	●
Patent Angle Indicator For Back-rest And Trendelenburg Position	●
φ113mm Wall Bumper Wheel	●
IV Pole	●
IV Pole Prevision	●
Drainage Hook	●
Three-grade Central Locking Pedal	●
ABS Mattress Retainer	●
Manual CPR	●
Electric CPR	●
Nurse Controller	●
Hand Controller	●
Weighing System	●
Night Light	●
Bed Extension 200mm	●
Trapeze Socket	●
Back Up Battery	●
Fifth and Sixth Siderails	●
Monkey Pole	○
Oxygen Cylinder Holder	○
Basket	○

Note: ●-YES, ○-Optional



CONSTANT INNOVATION

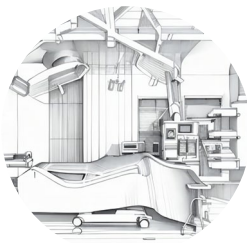


Design



Research



Manufacturing

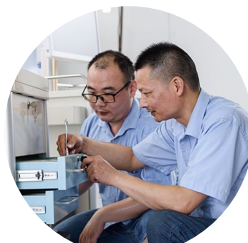

Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited
Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC
I7800 CASTLETON ST STE 180.
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue