

CELEBRATING
ANNIVERSARY

23

· 2002-2025 ·



Healthier your life

Z7z


ICU SMART BED

Registration code: SKD-B

www.saikangmedical.com

SAIKANG Medical®

THE EVOLUTION OF MEDICAL FURNITURE

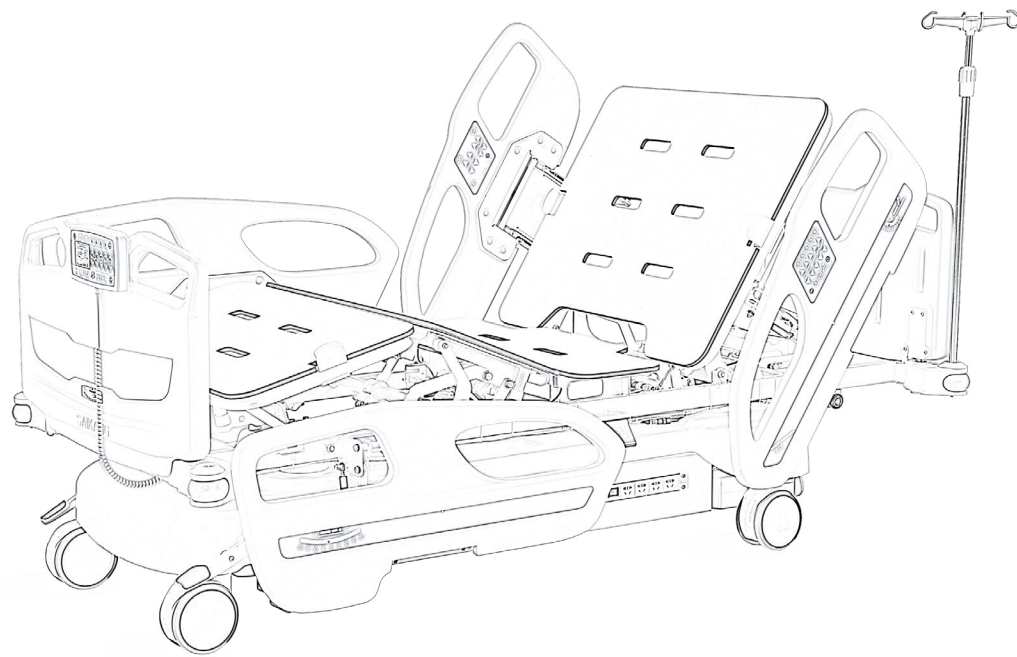

The unique germ-free design meets hospitals' strict infection control requirement. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization.

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

Panasonic Robotic ensure 360° full smooth welding.

RACING AGAINST TIME HEALTHIER YOUR LIFE

ENGINEERED FOR PATIENT WELLNESS, DESIGNED FOR CLINICAL EXCELLENCE

Working alongside individuals in critical situations, we recognized the pivotal role of the bed.

- The question arose: to what extent can care be enhanced by improving the operator's capabilities through the bed?

- Transforming the approach, we see designing a bed tailored to the caregiver as not a constraint but an influential and adaptable instrument.



“DESIGNED TO ASSIST ME IN
WORKING MORE EFFICIENTLY”



DUAL RETRACTION PATENTED DESIGN INNOVATIVE 75+15 mm



Auto-Regression

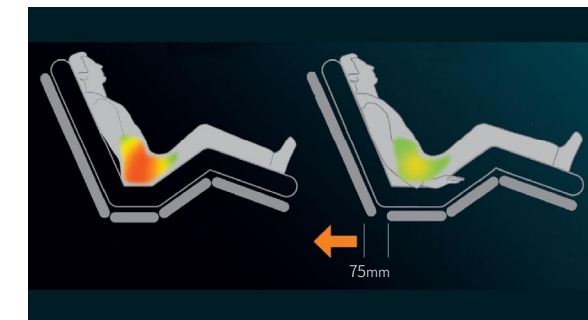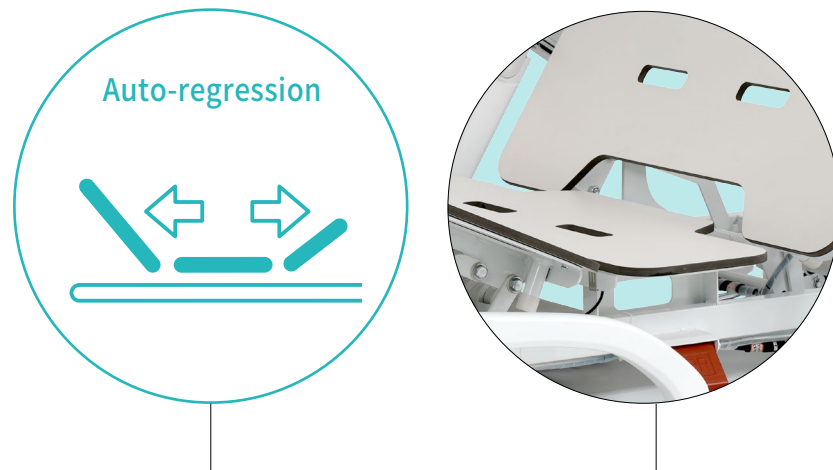

Pressure ulcers, or bedsores, result from prolonged pressure and friction on skin and underlying tissues.

- The dual retraction feature in a medical bed extends and distributes pressure evenly, reducing pressure concentration on the pelvic area.

- This design enhances patient comfort and prevents pressure ulcer occurrence proactively.

- The patented feature is applied to the backrest and thigh rest, retracting them by 75mm and 15mm respectively.

- Retracting these components extends the bed surface in the pelvic area, alleviating pressure and shear forces, thus minimizing pressure ulcer risk.

- The innovative design prioritizes patient well-being and offers a proactive solution for pressure ulcer prevention.

FULL-BODY X-RAY TRANSLUCENT BED PLATFORM



12mm thick phenolic resin bed platform for full-body X-ray functionality

Enlarge the seat area when back part is raised, providing more comforts for the patients.

Reduce the depression of the sponge at the seat panel.



Patient strap holes.

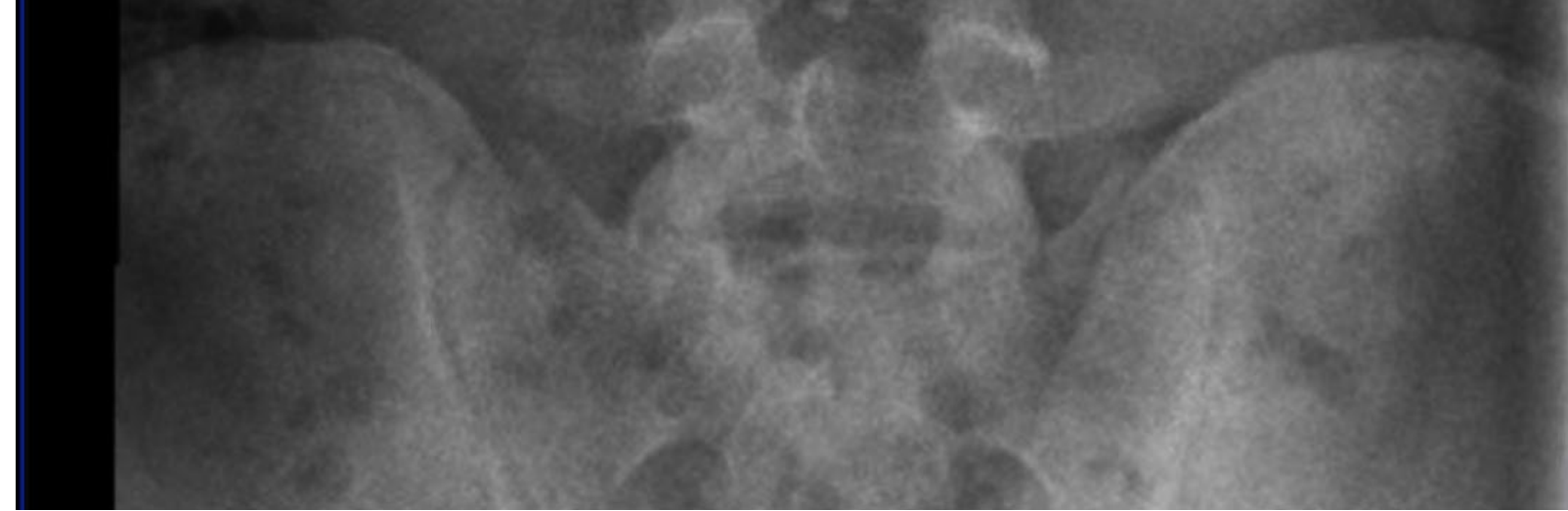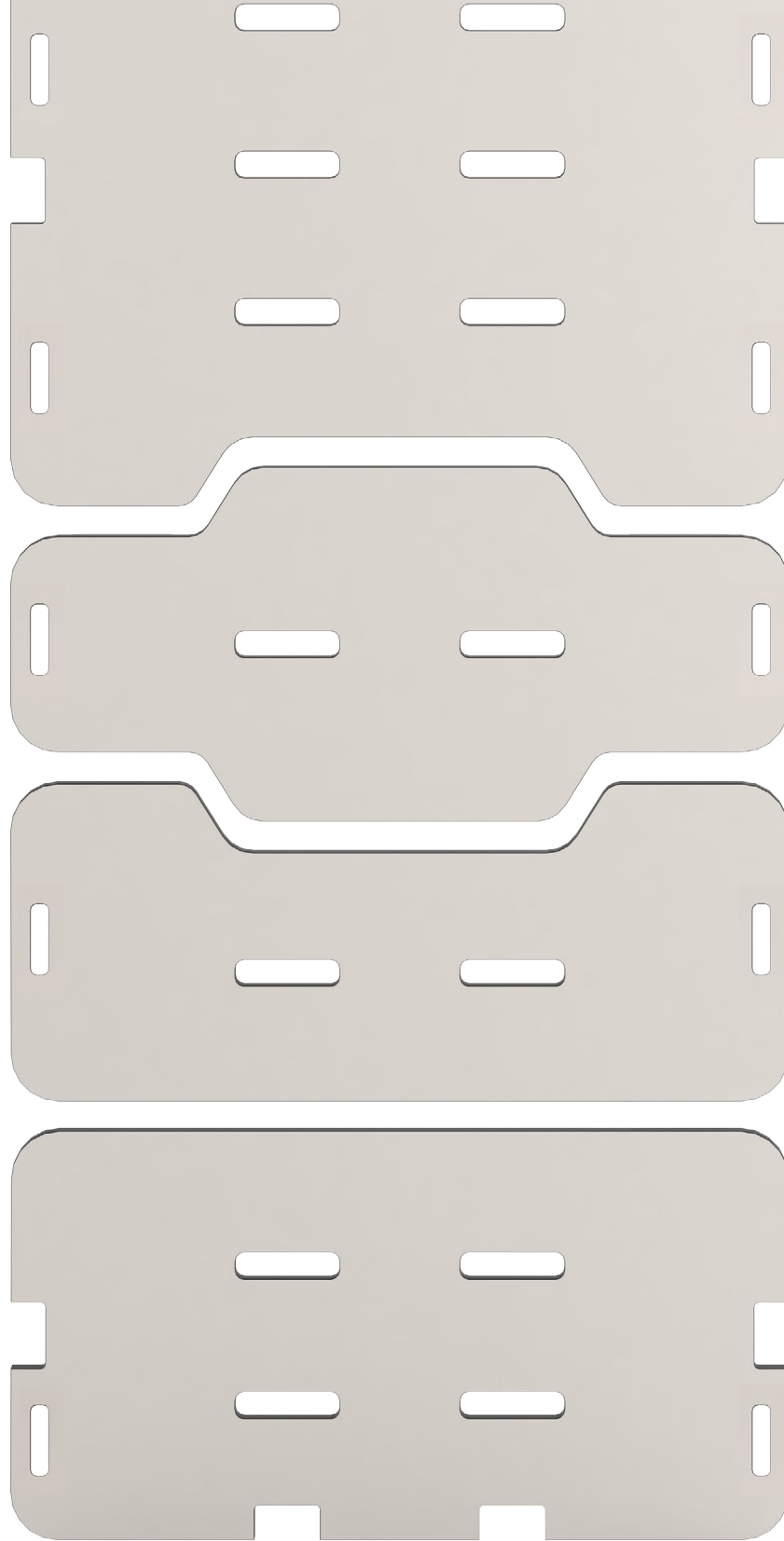

X-RAY EXAMINATION WITHOUT MOVING THE PATIENT


With 500*1800mm full-body X-ray translucence

Movable cassette holder size: 400*520mm

Wide area of radiotransparency.

It has an X-Ray cassette holder for full-body x-ray examinations, without moving the patient.

The X-ray tray holder fixed at the back section, easy for X-ray cassette film horizontally put out & put in, to avoid patient movement.

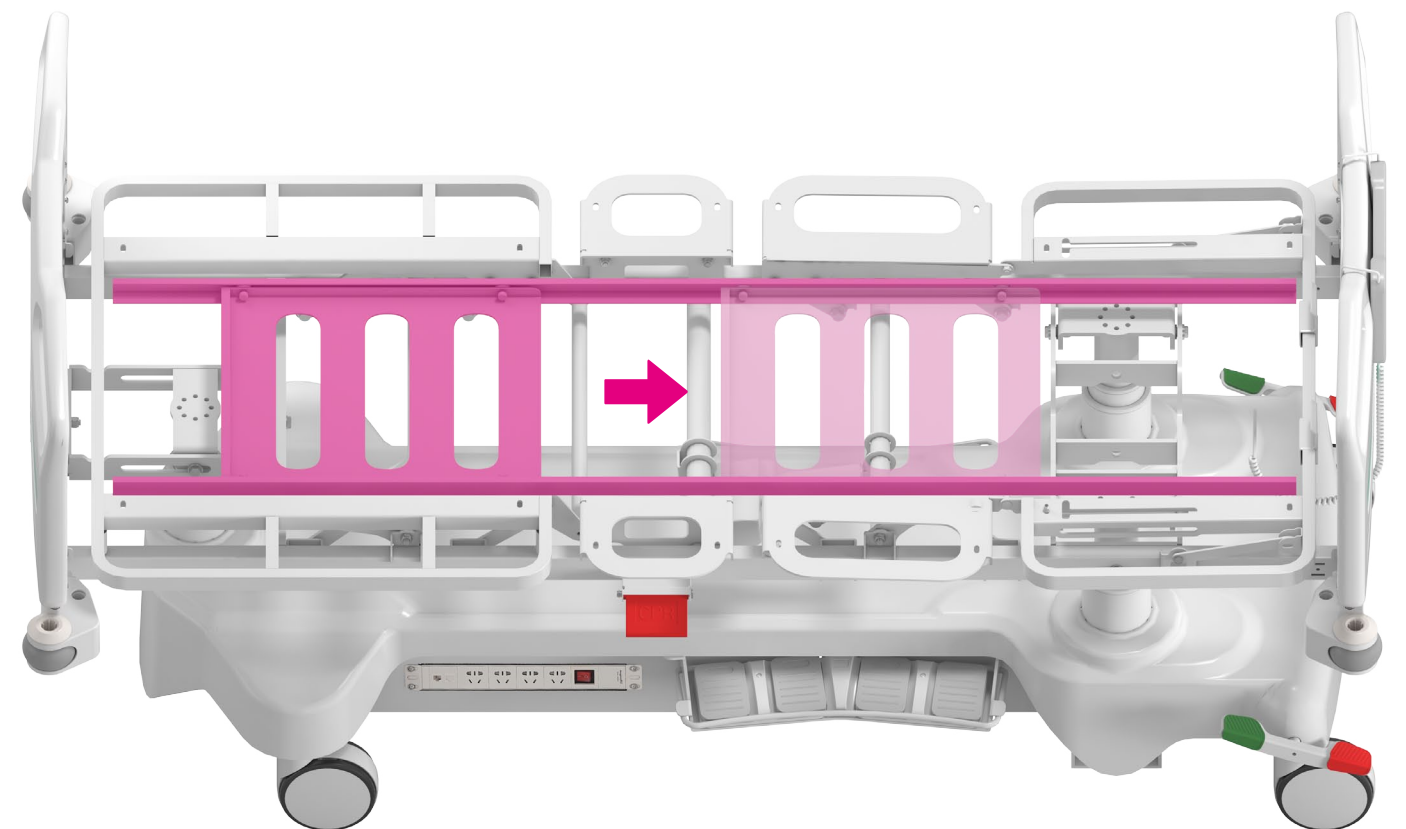



Back-rest



Bed height min.



Reverse Trendelenburg



Lateral tilt right



Knee-rest



Bed height max.



Trendelenburg



Lateral tilt left



The CPR adjustment allows to bring the mattress platform immediately down to the horizontal.



! FULL-COVERAGE SAFETY SIDE RAILS

Patient safety and protection are vital. Effective prevention can avoid fall injuries. When one rail is lowered, the bed's turning function locks to prevent turning, ensuring patient safety from falls due to unique side rail safety measures.

CUSTOMISED CONTROL CONCEPTS

Nursing staff and patients require safe and convenient bed operation options. Our 'light' version control panel provides customized features both inside and outside the bed. By prioritizing essential functions, operation becomes significantly simpler. An easily accessible handset with ergonomic buttons is also provided as an alternative.

Easy operation using the control panel “light”

Staff and patients benefit from the tailor-made convenience and care functions that are available on either side of the "light"-version control panel. The panel can be integrated in the split safety sides if required.




Glasses holder



Angle indicator



Backrest 30° indicator line

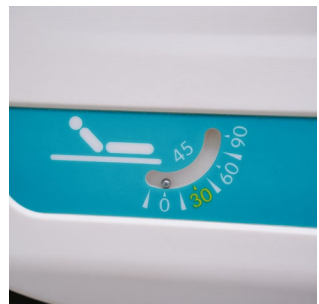

Backrest angle display at 30° and 45°

Just one-hand operation required to effortlessly release and lower the split safety siderail.

All it takes for the patient to get in and out of bed unhindered is to lower one safety siderail element in the centre of the bed.



Easy bed adjustment with the “control panel light” version

The “control panel light” version is always within easy reach of patients. The straightforward adjustment options are easy to understand even for people with cognitive impairment.

The “light”-version control panel is always at eye level for the patient.

PP BED ENDS WITH SAFE LOCK

Head and foot tilt angle
indicator

Integrated mechanical tilt angle
indicator on headboard and
footboard.



QUICK EMERGENCY ADJUSTMENTS

During resuscitation, it serves as a crucial aid. The shock position can be instantly activated using the easily accessible emergency lever at the head end or by using the control devices. Quickly removable headboards and footboards allow unimpeded access to the patient.



The orange emergency buttons on the outer side of the "light"-version control panel enable staff to act quickly without taking their eyes off the patient.



The CPR position can be set even faster using the mechanical emergency lever under the mattress base.




Effective use of bed space is crucial. The "integrated bumper wheels" design with IV pole holes and universal D series $\Phi 85\text{mm}$ bumper wheels provide effective collision protection during bed transport.



SS IV pole, bearing 15kg.



The headboard and footboard can be unlocked and removed without tools to allow improved access to the patient.

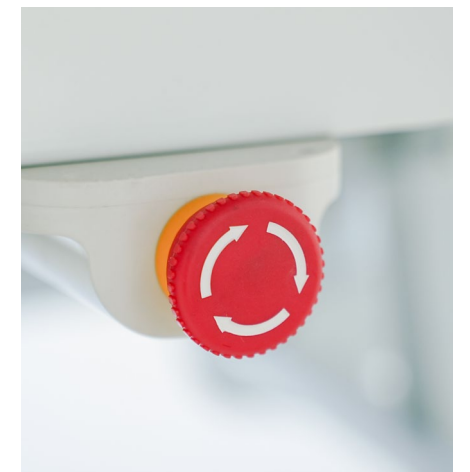

One button stop, prevent the maloperation.

LG Chem fresh PP material, plug in type.

FUNCTIONAL FEATURES




Each side has a set of foot pedals controlling bed height and tilts, with anti-mistake guard bar to prevent accidental activation.



6" Double side castors, central locking. TPR tire no worn out after running 10 km.
Anti-winding hard shell cover, united forming without bolts.
Pass dynamic test: pass obstacles 40000 times, under the standard 200kg load of the whole bed.



One multi-outlet power strip (4 outlets) on each side, providing power for various ICU instruments.

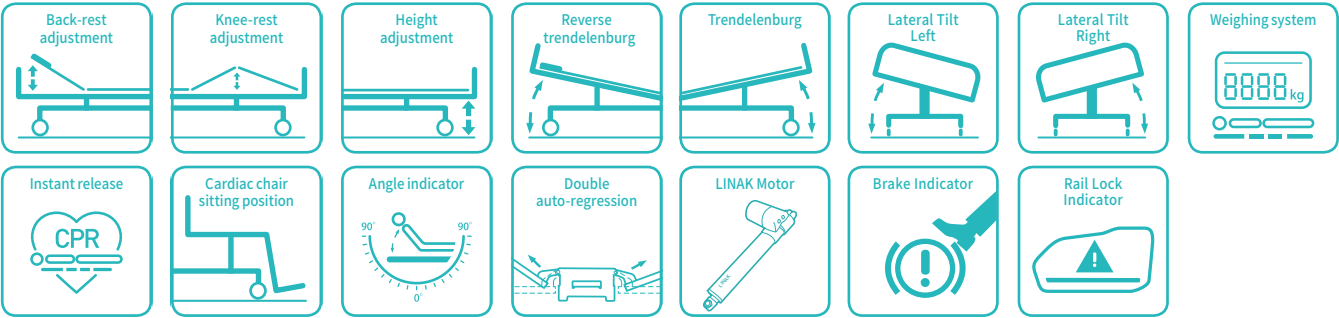


PRODUCT INFORMATION

Model	Z7z
External size (LxW)	2127x1070 mm
Mattress platform	1925x900 mm
Safe working load	250 kg
Height adjustment	530-930 mm
Back-rest adjustment	0-75° (±10°)
Knee-rest adjustment	0-35° (±10°)
Reverse trendelenburg	0-15° (±2°)
Trendelenburg	0-15° (±2°)
Lateral Tilt	0-17° (±2°)
Angle indicator	0-90°
Auto-regression	75+15 mm
Position line	30°
LINAK motor	●
Telescopic column motor	●
6" Double side secure® castors, central locking	●
PP side rall with integrated LCD control panel	●
PP bed ends with safe lock	●
X-ray transparent phenolic resin bed platform	●
Patent angle indicator	●
Φ85mm wall bumper wheel	●
Trapeze socket	●
IV pole	●
IV pole prevision	●
Drainage hook	●
Central locking pedal	●
PP mattress retainer	●
Nurse Cotroller	●
Emergency stop button	●
Night light	●
Manual CPR	●
Electric CPR	●
Weighing system	●
Monkey pole	○

Note: ●-YES, ○-Optional



CONSTANT INNOVATION

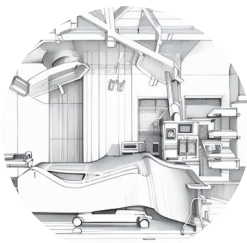


Design



Research



Manufacturing

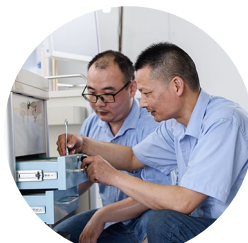

Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited
Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC
17800 CASTLETON ST STE 180,
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue