



A307

ELECTRIC-HYDRAULIC OPERATING TABLE

Registration code: SKT-C



SAIKANG Medical

Healthier
your life

A307

FEATURES AND HIGHLIGHTS

Technical parameters

▪ Length/Width:	2140 × 550 mm ± 10 mm
▪ Height Elevating/Lowering:	550-1050 mm ± 10 mm
▪ Forward/Backward Tilt:	30° ± 3°
▪ Left/Right Tilt:	20° ± 3°
▪ Backrest (Up/Down):	80° ± 3° / 50° ± 3°
▪ Head Section (Up/Down):	45° ± 3° / 90° ± 3°
▪ Leg Section (Up/Down):	15° ± 3° / 90° ± 3°
▪ Leg Section Outward:	90° ± 3°
▪ Horizontal Sliding:	300 mm ± 5 mm
▪ One-Button Flexion:	220° ± 3°
▪ Reverse Flexion:	100° ± 3°

Function Configurations

▪ Electric Brake
▪ Dual Head Section
▪ Dual Articulated Leg Section
▪ Head/Leg Section Interchangeable with Automatic Recognition

Technical configuration

☑ Anesthesia Screen	1 unit
☑ Shoulder Supports	1 pair
☑ Arm Rests	1 pair
☑ Lumbar Supports	1 pair
☑ Leg Supports	1 pair

Design

This electric hydraulic operating table is suitable for thoracic and abdominal surgery, brain surgery, ENT surgery, obstetrics and gynecology, urology, orthopedics and other surgical operations.

The material is 304 stainless steel, and the platform is made of high-strength synthetic plate which is translucent for X-ray. It can also shift the table back and forth up to 300mm for C-arm examination of all parts of the patient's body during the operation without moving the patient.

Welding | Craftsmanship

Adopted robot integral welding, smooth appearance, stable and reliable.



Material

The whole machine is made of 304 stainless steel material, features beautiful, strong, impact resistant, damage resistant, corrosion resistant and more lightweight. And with strong bearing capacity of 360kg.

Base

Base is made of Q235 board, with 304 stainless steel cover which is anti-corrosion and easy to clean.

Surface

Five-section modular table design, the head plate, upper back plate, and leg plate can be disassembled manually. it can meet different surgical needs by increasing or reducing the table modules.

Mattress

Memory sponge mattress (polyurethane memory foam), thickness 75mm, seamless mattress surface, waterproof, air permeable, easy to clean, anti-static, and can be see-through.

Head plate

The head plate is a double-jointed head plate that allows the patient to obtain a comfortable head position during the surgery.

Back plate

Divided into upper backplate and lower backplate, the upper backplate can be disassembled by manually, easy to install.

Leg plate

The leg plate is a double-joint horizontal rotating leg plate, so there's better operating space for surgeon.

Drive mode

Adopted electro-hydraulic control mode, flexibly control all table movements through the micro-key switch controller, and electrically adjust 7 main action groups including bed lifting, reversed trendelenburg and trendelenburg, left and right tilt, backplate lifting, leg plate lifting, countertop translation, and brake.

Hydraulic control system

Adopted domestic well-known hydraulic control system (hydraulic pump, solenoid valve, hydraulic cylinder), with strong reliability, good stability, no frustration, quiet and smooth.

Modular tabletop

The head plate and leg plate can be interchanged according to the operation requirements, and the hook type one-key quick disassembly structure is adopted. The system can automatically identify the direction of the head and leg plate after interchange.

Horizontal sliding

The table has a translation function which is driven by an independent motor, and the translation distance is $\geq 300\text{mm}$.

One button reset

The one-button reset button can complete the horizontal, forward and reverse flexion functions with just one key.


Controller

Functions controlled by wired hand controller or column control panel. Each control mode is independent, and the dual-mode operation design provides double security for the operation.



Buckling & reverse buckling

With functions of one-key buckling $\geq 220^\circ$ and one-key reverse buckling $\leq 100^\circ$.

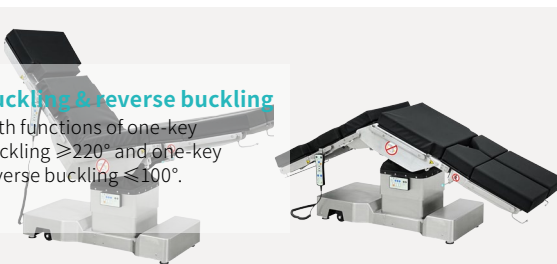

Emergency stop button

The emergency stop button is added to further improve the safety of operating bed.




Accessory Quick Disconnect Device

Accessories adopt a structure of quick-disassemble device, flexible and convenient.



Castor

Using industrial castors, moving quiet, smooth and convenient. Diameter 75mm. Four casters are all universal wheels.

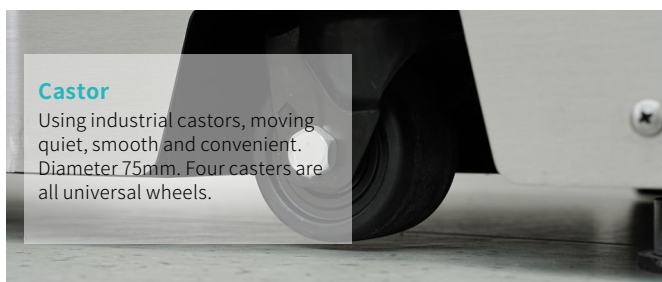


Battery

Equipped with high-performance rechargeable batteries, it can meet the needs of ≥ 50 surgeries, ensuring that the operating bed can work without AC power supply, ensuring maximum safety and reliability. Rechargeable batteries require no maintenance and can be used for a long time.




Anesthesia screen frame

Equipped with a stainless steel anesthesia screen frame for hanging anesthesia curtains, used to shield and cover the patient's head during use.

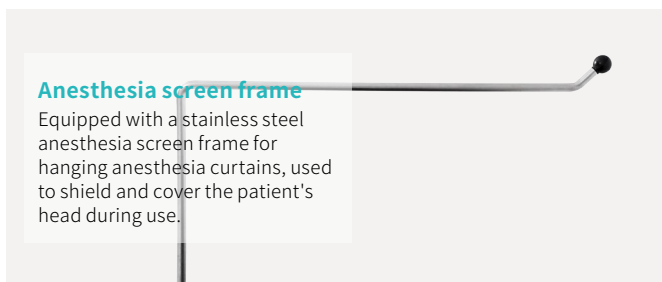


Shoulder support

Polyurethane support frame, adjustable in angle, used to stabilize the patient's body.

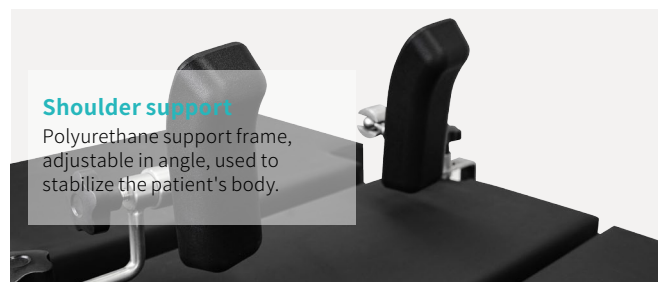


Body support

Polyurethane support frame, adjustable in angle, used to stabilize the patient's body.

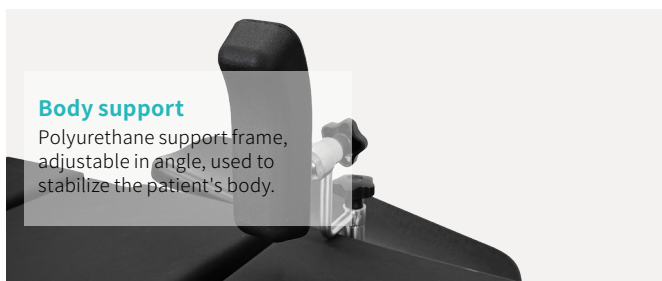


Arm rest

Polyurethane widened armrest board with attached straps, providing stable performance and versatile applications.

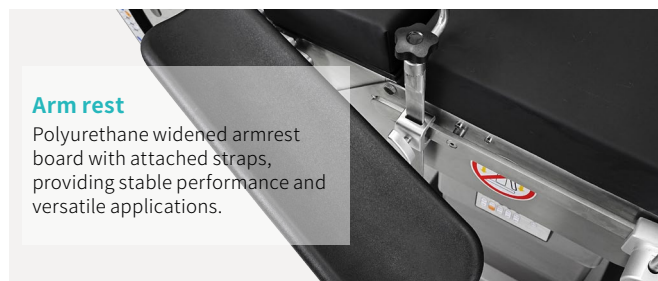


Leg rest

Manual multi-angle adjustment for easy use.

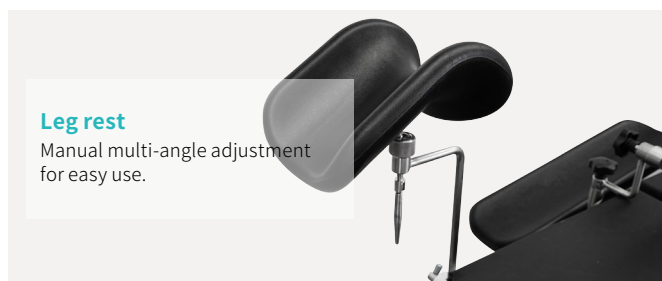

Foot plate

Be used to support patient's feet.

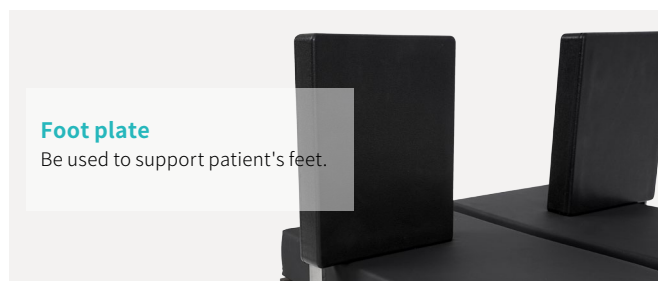

Professional functions

Lowest position of table top
550mm



Highest position of table top
1050mm



Back plate upward
0°-80°



Back plate downward
0°-50°



Reversed trendelenburg
0°-30°



Trendelenburg
0°-30°



Left lateral tilt
0°-20°



Right lateral tilt
0°-20°



Head plate up
0°-45°



Head plate down
0°-90°



Leg plate up/down
15°/90°



Leg plate outwards
0°-90°




Horizontal sliding
0-300mm





CONSTANT INNOVATION

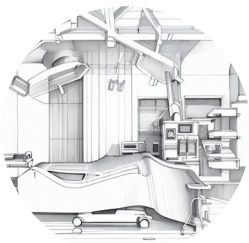


Design



Research



Manufacturing

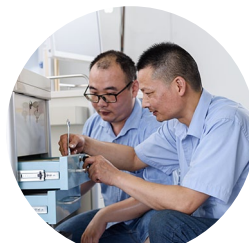

Maintenance



Service



WORLDWIDE BRANCHES

Portugal Showroom

Av. da Liberdade, nº 101
8150-101
São Brás de Alportel
Portugal
TEL: +351 289 090 212

Brazil Showroom

Rua Lauro Muller 116, Sala
1503 - Botafogo, Rio de
Janeiro - RJ, 22290-160
E-mail: mehecobrasil@
gmail.com

Indian Showroom

1460 M.I.E, Bahadurgarh,
Haryana-124507, INDIA
TEL: +91 96677 70747
E-mail: ind@saikangmedical.
com

China Headquarter

No.35 Lehong Road,
Zhangjiagang City,
Jiangsu Province, China.
TEL: +86 18021231901
E-mail: export@
saikangmedical.com

South East Asia Branch

**SAIKANG Medical
(Thailand) Co., Ltd.**
Ample Tower 909 Room 11/3
Debaratana Road,
Bangna Nuea, Bangkok
10260 Thailand
TEL: +66-936-569469

US Branch

**Healthward Technology,
INC**
17800 castleton st ste 180,
Los Angeles
TEL: +1 (380) 254 4913

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue