

CELEBRATING
ANNIVERSARY

23rd

· 2002-2025 ·



Healthier your life

V8v

ELECTRIC ICU BED

Registration code: SKD-CV

www.saikangmedical.com

SAIKANG Medical®


THE EVOLUTION OF MEDICAL FURNITURE

The unique germ-free design meets hospitals' strict infection control requirement. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization.

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

Panasonic Robotic ensure 360° full smooth welding.

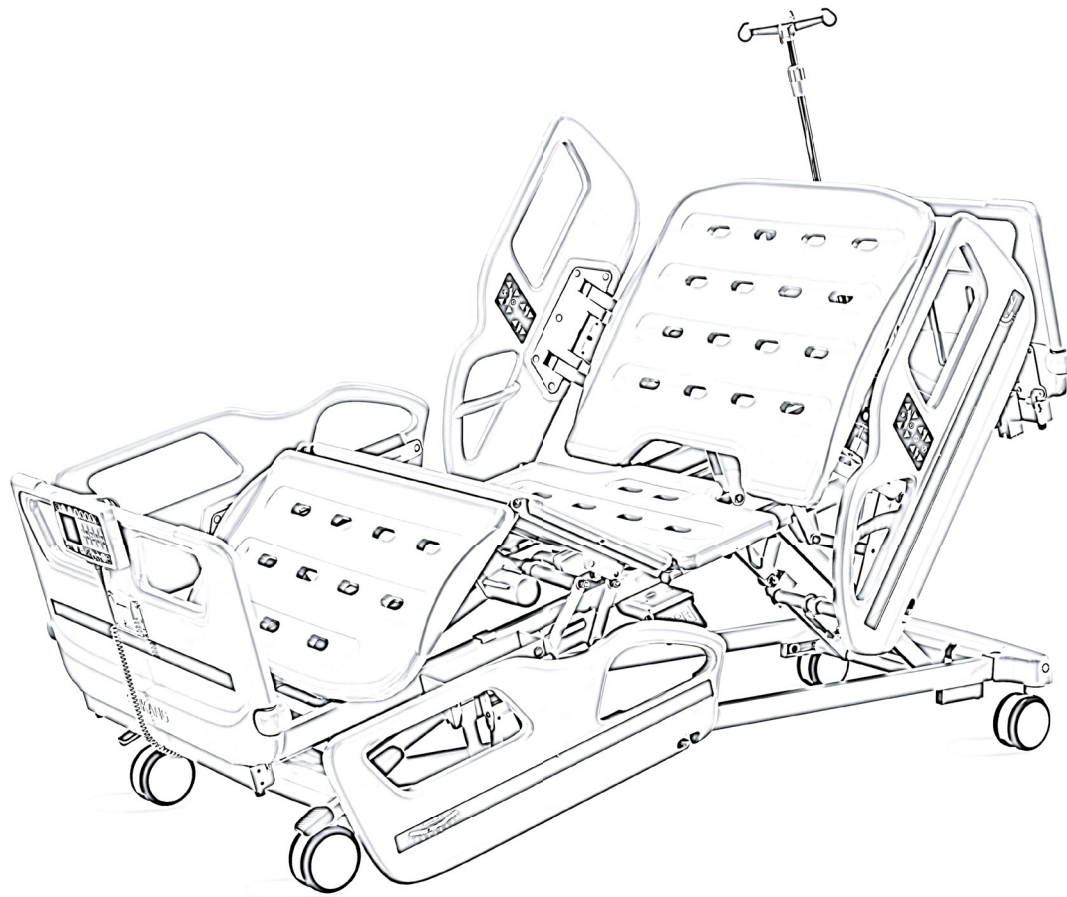

ENGINEERED FOR PATIENT WELLNESS, DESIGNED FOR CLINICAL EXCELLENCE

Working alongside individuals in critical situations, we recognized the pivotal role of the bed.

- The question arose: to what extent can care be enhanced by improving the operator's capabilities through the bed?

- Transforming the approach, we see designing a bed tailored to the caregiver as not a constraint but an influential and adaptable instrument.

“DESIGNED TO ASSIST ME IN
WORKING MORE EFFICIENTLY”





RACING AGAINST TIME HEALTHIER YOUR LIFE




Backrest 30-degree smart stop: The clinically recommended position for preventing respiratory complications. Research shows that the 30-degree position significantly reduces the risk of lung infections.

Auto-Regression

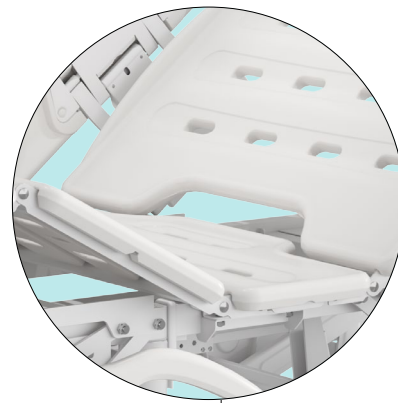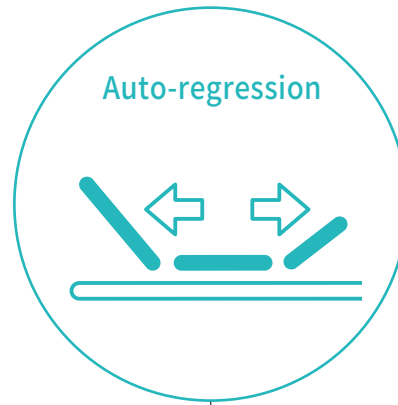

DUAL RETRACTION PATENTED DESIGN INNOVATIVE 67+47 mm

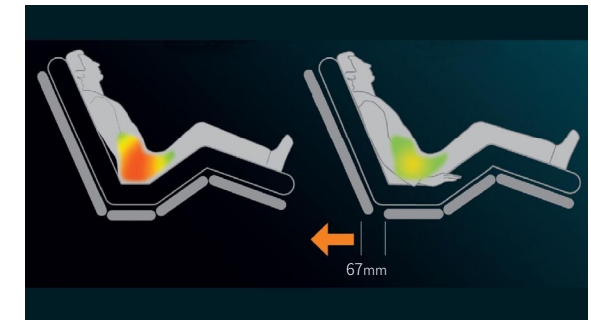

Pressure ulcers, or bedsores, result from prolonged pressure and friction on skin and underlying tissues.

- The dual retraction feature in a medical bed extends and distributes pressure evenly, reducing pressure concentration on the pelvic area.

- This design enhances patient comfort and prevents pressure ulcer occurrence proactively.

- The patented feature is applied to the backrest and thigh rest, retracting them by 67mm and 47mm respectively.

- Retracting these components extends the bed surface in the pelvic area, alleviating pressure and shear forces, thus minimizing pressure ulcer risk.

- The innovative design prioritizes patient well-being and offers a proactive solution for pressure ulcer prevention.

INTEGRATED PP BED PLATFORM

Integrated PP bed platform, adopting advanced manufacturing technology.


4-5 section platform, enables independent installation and easy maintenance.

Inlarge the area of the seat panel when back part is raised, providing more comforts for the patients.

The seat compensation part reduces sponge depression on the seat panel.




The bed is equipped with mattress retainer all around the board, so that the mattress will not easily slide off.

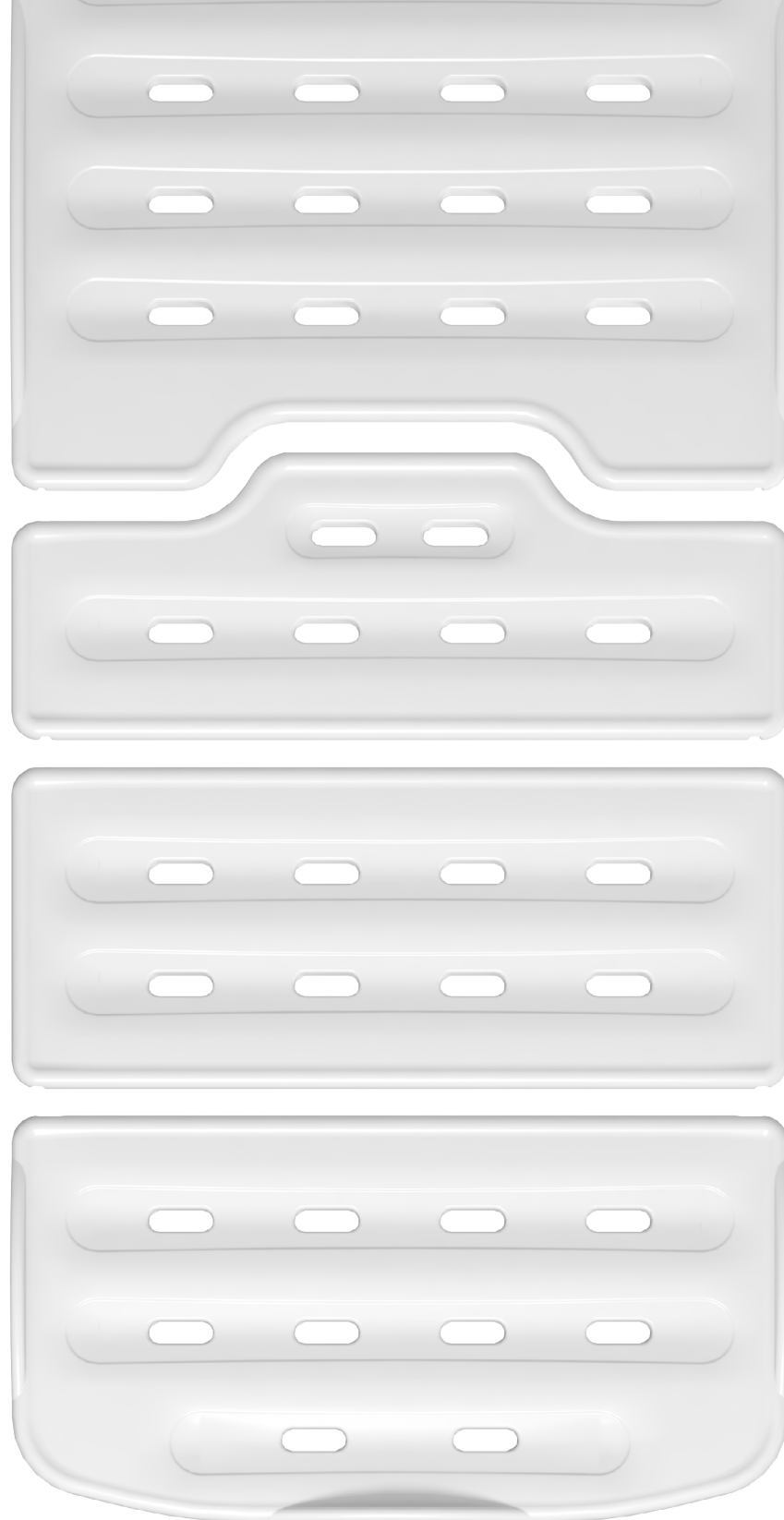

X-RAY EXAMINATION WITHOUT MOVING THE PATIENT

Wide area of radiotransparency. The PP radiotransparent backrest has an X-Ray cassette holder for chest x-ray examinations, without moving the patient and at any angle of the backrest.

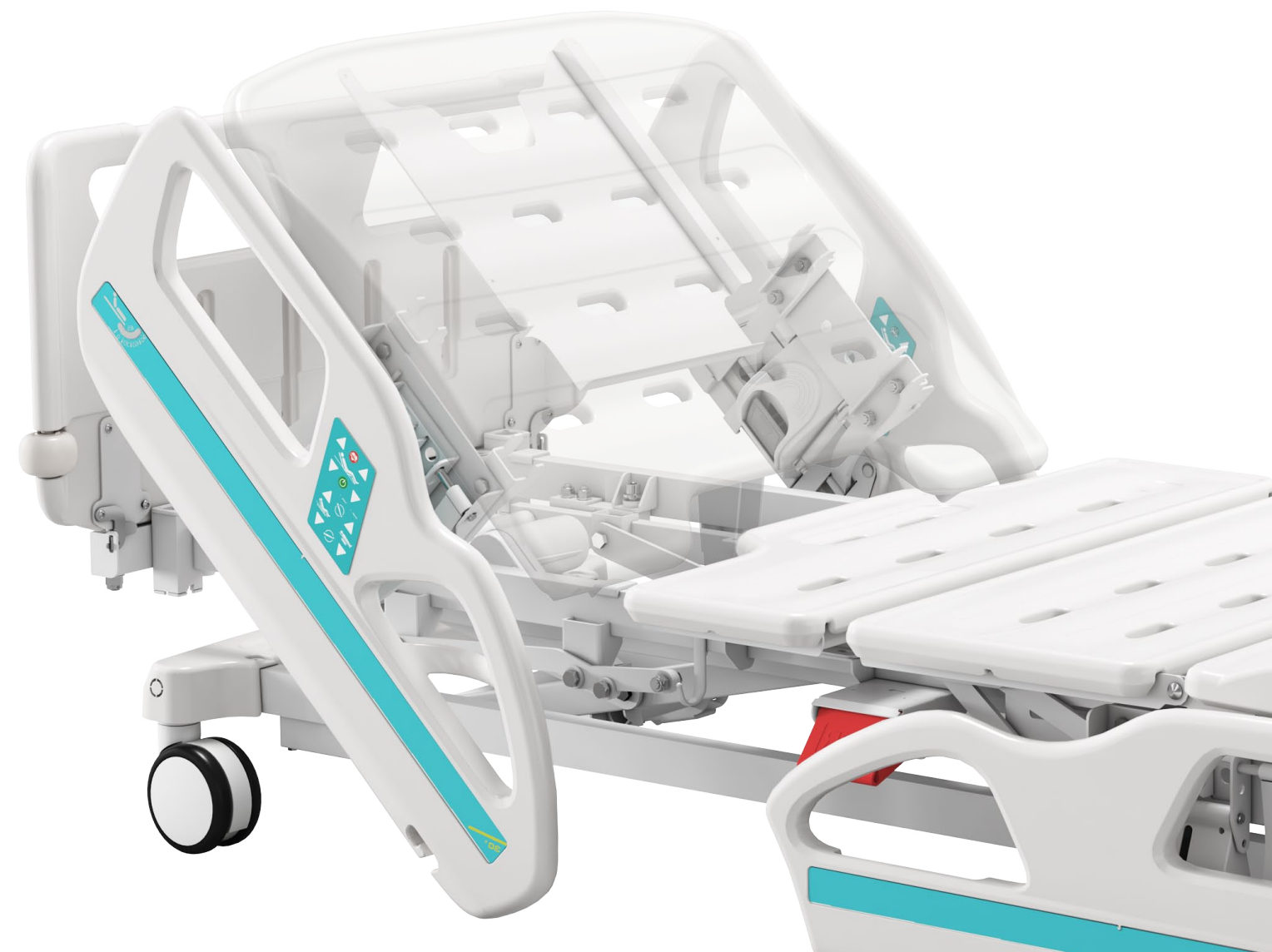



① Back-rest



③ Bed height min.



⑤ Reverse Trendelenburg



⑦ The examination position height adjustment facilities medical staff in the examination, medication and monitoring operations, by standing up straight.



② Knee-rest



④ Bed height max.



⑥ Trendelenburg



⑧ The CPR adjustment allows to bring the mattress platform immediately down to the horizontal.



⑨ The chair position helps prevent respiratory complications, it induces to a more physiologic position and has a psychologically positive influence.

ADVANTAGES FOR NURSING CARE

Tuck away side rail with buffer.

Built-in controller in the head guardrail, allowing the patient to adjust the basic position of the bed easily.

Integrated handle facilitates the patients to get up/off the bed.

SAFE AND EASY TO OPERATE



CUSTOMISED CONTROL CONCEPTS

Nursing staff and patients require safe and convenient bed operation options. Our 'light' version control panel provides customized features both inside and outside the bed. By prioritizing essential functions, operation becomes significantly simpler. An easily accessible handset with ergonomic buttons is also provided as an alternative.

Easy operation using the control panel “light”

Staff and patients benefit from the tailor-made convenience and care functions that are available on either side of the "light"-version control panel. The panel can be integrated in the split safety sides if required.



Just one-hand operation required to effortlessly release and lower the split safety siderail.

All it takes for the patient to get in and out of bed unhindered is to lower one safety siderail element in the centre of the bed.



SAIKANG Medical®

Easy bed adjustment with the “control panel light” version

The “control panel light” version is always within easy reach of patients. The straightforward adjustment options are easy to understand even for people with cognitive impairment.

The “light”-version control panel is always at eye level for the patient.



OPERATION WITH ERGONOMIC HANDSET

The function keys with a dark blue background are easy to understand in any situation.



The intuitive handset is available for both types of safety side and can always be attached conveniently and within easy reach.



The ergonomic handset with its intuitive symbols seems familiar to every patient straight away and creates a positive feeling of reassurance.



QUICK EMERGENCY ADJUSTMENTS

During resuscitation, it serves as a crucial aid. The shock position can be instantly activated using the easily accessible emergency lever at the head end or by using the control devices. Quickly removable headboards and footboards allow unimpeded access to the patient.



The orange emergency buttons on the outer side of the "light"-version control panel enable staff to act quickly without taking their eyes off the patient.



The CPR position can be set even faster using the mechanical emergency lever under the mattress base.



The headboard and footboard can be unlocked and removed without tools to allow improved access to the patient.

LG Chem fresh PP material, plug in type.


Equipped with the bumpers wheel, diameter 65mm, preventing collision.



One button stop, prevent the maloperation.

EFFORTLESS MOBILITY: SMOOTH AND SECURE REPOSITIONING WITH MINIMAL CAREGIVER STRAIN

Its compact dimensions make it ideal for use in multi-bed rooms. Good maneuverability is crucial to facilitate the work of staff and prevent collisions. The modern castors make the bed easy to move, while an optional fifth castor facilitates maneuvering around tight corners.

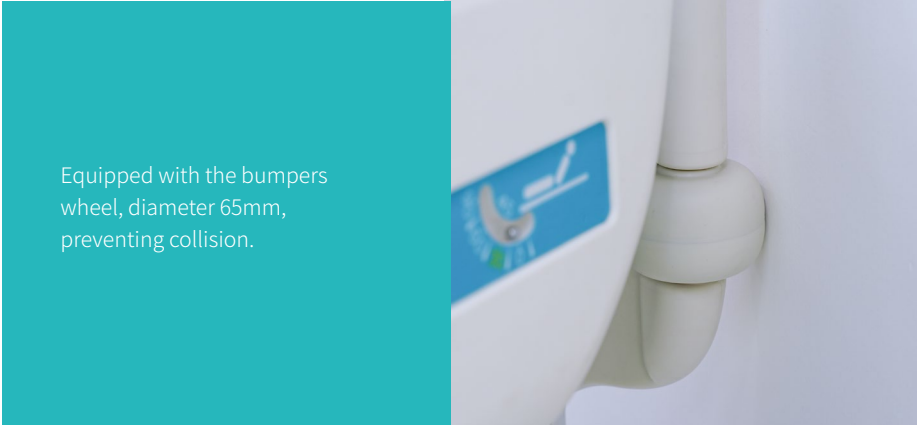

Equipped with the bumpers wheel, diameter 65mm, preventing collision.



6-inch central locking double-sided castors. TPR tire no worn out after running 10 km.
Anti-winding hard shell cover,united forming without bolts.
Pass dynamic test: Successfully pass 40,000 obstacle tests without deformation, under the standard 200kg load of the entire bed.
1-speed unlock; 1-speed universal, 1-speed directional.



SAFE HYGIENE

Easy manual reprocessing

During a pandemic, safe and verifiable hygiene is more crucial than ever. Its low-recess design simplifies manual reprocessing for the hygiene team, saving time, physical effort, and personnel resources.



The lowest height position minimises the rate of fall in the event of any fall. This significantly reduces the risk of injury.



The split safety side element at the head end acts as a mobilisation support to help the patient get out of bed.




No inaccessible recesses in the headboard or footboard and safety sides..



SS IV pole, bearing 15kg.



Plastic drainage hooks, movable.

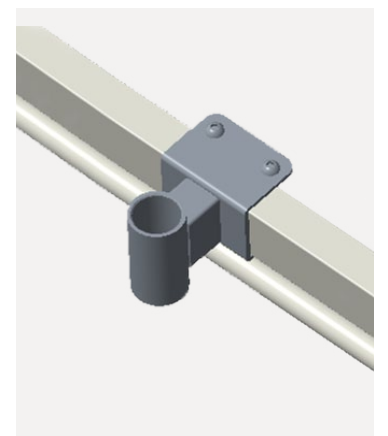

Trapeze Bracket



Level mattress base with easy to remove plastic cover.



Tidy cable routing.



The nurse control panel features a flexible spiral cable, allowing easy and flexible movement to operate the bed in the optimal position. The nurse can control all standard functions and various one-touch automatic positions, including the cardiac chair position.



The under bed light helps patients by lighting their way around the room at night.



Back up battery, 18 times movement with 100kg bearing, 130 times movement without bearing.

PATIENT WEIGHT UNDER CONTROL AT ALL TIMES




The exclusive weighing system was designed to guarantee a high level of precision independently from the mattress platform's position.


Thanks to the 4 sensors located at the corners of thebed, the load is better distributed and the reading level more accurate.

The weighing data can be recorded manually whenever necessary or automatically every 24 hours.


While the weight of the patient is monitored it is possible to set up different types of alarm: exceededthreshold weight alarm or “bed exit” alarm.

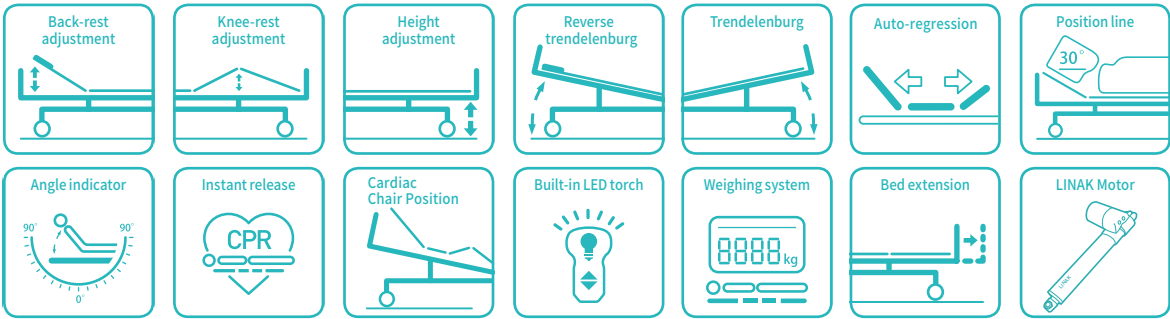

Options

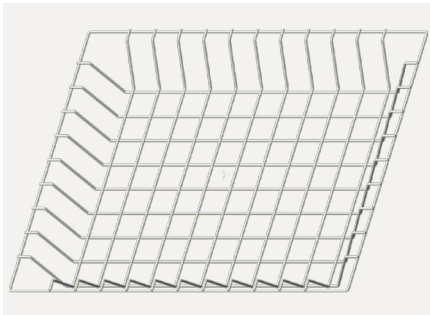


Utility Shelf
Detachable unitily shelf for patients' personal stuff.

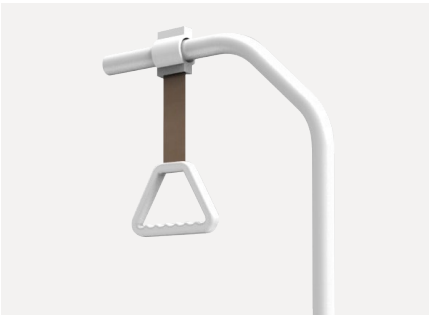

Trapeze Pole
32mm dia heavy duty trapeze pole.



Mechanical Lever

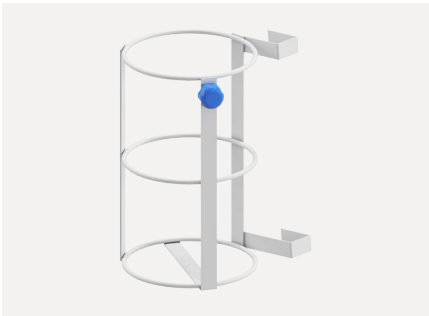

Oxygen Cylinder Holder

Optional color (MOQ100pcs- color changing)




Product Information

Model	V8v
External size (LxW)	2085x1025mm±10mm
Mattress platform	1925x900mm±10mm
Safe working load	250kg
Height adjustment	420-870mm±10mm
Back-rest adjustment	0-75° (±10°)
Knee-rest adjustment	0-28° (±10°)
Reverse trendelenburg	0-17° (±2°)
Trendelenburg	0-17° (±2°)
Angle indicator	0-90°
Auto-regression	67+47mm
Position line	30°
LINAK Motor	●
6" Double Side Secure®Castors, Central Locking	●
PP Platform with Integrated Mattress Retainer	●
PP Side Rail with Integrated Control Panel	●
PP Bed Ends plug-in type, with safe lock	●
Patent Angle Indicator for backrest and Trendelenburg Position	●
IV Pole 4 hooks	●
IV Pole Prevision	●
Drainage Hook	●
Central Locking Pedal	●
Hand Cotroller	●
Nurse Cotroller	●
Back-up Battery	●
Cardiac Chair Position	●
Electric CPR	●
Manual CPR	●
Auto Break at 30°	●
Emergency Stop Button	●
Night Light	●
Cassette Holder (Backrest)	●
Bed Extension	●
Trapeze Socket	●
Weighing system	●
Monkey Pole	○
Oxygen Cylinder Holder	○
Basket	○
Mechanical Lever	○

Note: ●-YES, ○-Optional



CONSTANT INNOVATION

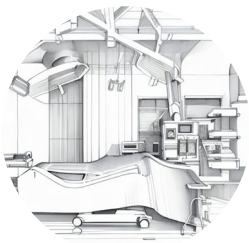


Design



Research



Manufacturing

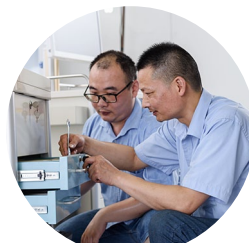

Maintenance



Service



WORLDWIDE BRANCHES

Portugal Showroom

Av. da Liberdade, nº 101
8150-101
São Brás de Alportel
Portugal
TEL: +351 289 090 212

Brazil Showroom

Rua Lauro Muller 116, Sala
1503 - Botafogo, Rio de
Janeiro - RJ, 22290-160
E-mail: mehecobrasil@
gmail.com

Indian Showroom

1460 M.I.E, Bahadurgarh,
Haryana-124507, INDIA
TEL: +91 96677 70747
E-mail: ind@saikangmedical.
com

China Headquarter

No.35 Lehong Road,
Zhangjiagang City,
Jiangsu Province, China.
TEL: +86 18021231901
E-mail: export@
saikangmedical.com

South East Asia Branch

**SAIKANG Medical
(Thailand) Co., Ltd.**
Ample Tower 909 Room 11/3
Debaratana Road,
Bangna Nuea, Bangkok
10260 Thailand
TEL: +66-936-569469

US Branch

**Healthward Technology,
INC**
17800 castleton st ste 180,
Los Angeles
TEL: +1 (380) 254 4913

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue