



Oasis Series

Y8t

MULTIFUNCTIONAL ELECTRIC BED

Registration code: SKD-C




SAIKANG/Medical

Healthier
your life

Y8t

FEATURES AND HIGHLIGHTS

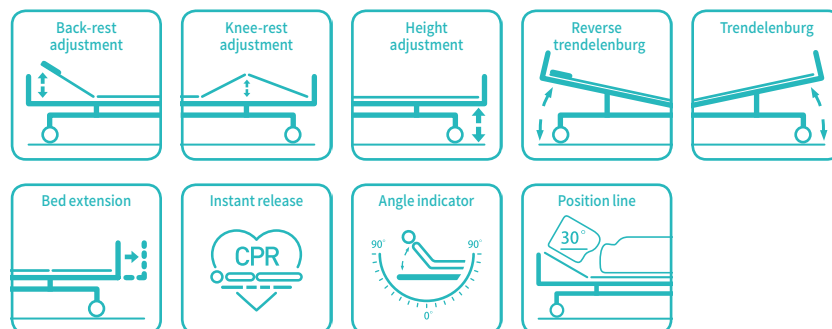

Technical parameters

▪ External size (LxW):	2140x1060 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	480-770 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Reverse trendelenburg:	0-12° (±2°)
▪ Trendelenburg:	0-12° (±2°)
▪ Angle indicator:	0-90°
▪ Position line:	30°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	

Optional

☐ Hand controller
☐ Weighing system
☐ Nurse control panel
☐ Cardiac chair position
☐ Bed extension 20cm
☐ Back-up battery
☐ Directional wheels function
☐ ABS mattress stopper
☐ Monkey pole
☐ Balkan frame
☐ Linen holder

Technical configuration

☒ LINAK motor	4 pcs
☒ 5" double side castors, central locking system	4 pcs
☒ PP Oasis side rail extra height	1 set
☒ PP Oasis bed ends with safe lock	1 set
☒ Stamping forming bed platform	1 pc
☒ Plastic socket inside pp guardrail for phone or glasses	2 pcs
☒ Patent angle indicator for back-rest and trendelenburg position	4 pcs
☒ Ø113mm wall bumper wheel	4 pcs
☒ IV pole	1 set
☒ IV pole prevision	4 pcs
☒ Drainage hook	2 pcs
☒ Central locking pedal	2 pcs
☒ Emergency stop button	1 pc
☒ Manual CPR	
☒ Electric CPR	

Go button

Use the GO button to activate the bed; protects the patient against an injury caused by any unintentional movement of the bed. Control system is automatically locked after 3 min.




One button mobilisation position

Respects easy patient mobilisation: One-button for patient mobilisation position -back plate raised, double hole design on back part guardrail, which supports patient self standing.




Lowest position indicator

The lowest position indicator design allows the medical staff to observe the bed safe position during night or long distance; green light indicate lowest bed height.

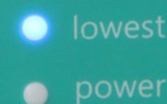

Great multi-care

Respects easy patient mobilisation:
One-button for patient mobilisation position -back plate raised, double hole design on back part guardrail, which supports patient self standing.

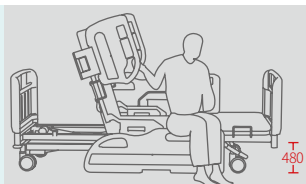**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding.
Solid bed base made of 5.0 mm stamped rotating plate.

**Steels | Painting**

Baosteel from Fortune Global 500.
11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Electric motor**

LINAK motor, with UL/ROHS/EN approval. Noise under 50dB, with lifetime 20,000times movement.

**Safety**

IEC60601 Safety Standard.
One-button emergency stop switch can quickly stop all going movement.

**CPR**

CPR emergency reset handle to immediately flatten the bed back part; located on both sides of the bed. Electric CPR button is setted inside control panel, one button to make the back plate flat.

**Castor**

5" diameter central locking castors with TPR tire, no worn-out after running 30KM, united forming without bolts.
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**Brake Padel**

The double-folded large-area rubber design, horizontal and vertical two-sided rubber, allows the nurse to touch the brakes maximally, saving effort.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Angle indicator**

Patent angle indicator, clearly show the inclination angles.

**Head & foot board**

Patented design Detachable head & foot board, using the eccentric locking structure, one step quick removable.

**Siderail**

4 pcs PP siderail (T20), with patent angle indicators on each one, clearly show back-rest and tredelenburg position.

**Siderail controller**

Control panel in the siderails:
Numerous studies have shown that the built-in controller located in the head guardrail is easier to reach, allowing the patient to adjust the basic position of the bed through the controller.

**Wall bumper wheel**

Equipped with the bumpers wheel, diameter 113mm.
Corner bumpers at the outermost periphery of the bed to protect the bed and wall

**IV pole**

The I.V. pole is made of stainless steel and has four plastic hooks that can bear up to 15 kg of weight.

**Drainage hooks**

2 pieces of plastic drainage hook.



Hand controller

Hand controller, keep qualified after 2000 times use.



Weighing system.



Night light under the bed





Nurse control panel, hook type.



One-button for patient mobilisation position -back plate raised, double hole design on back part guardrail, which supports patient self standing.




One button nursing position, helping the nurse to give fast treatment in a very comfortable height.

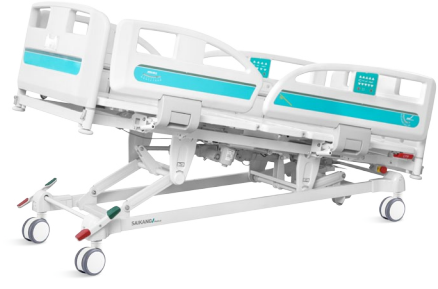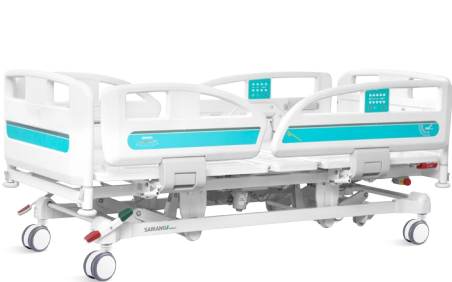



CONSTANT INNOVATION

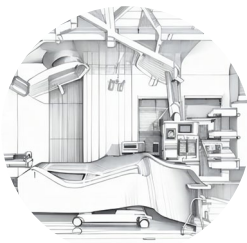


Design



Research



Manufacturing

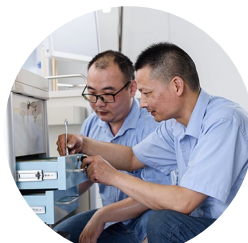

Maintenance



Service



WORLDWIDE BRANCHES

Portugal Showroom

Av. da Liberdade, nº 101
8150-101
São Brás de Alportel
Portugal
TEL: +351 289 090 212

Brazil Showroom

Rua Lauro Muller 116, Sala
1503 - Botafogo, Rio de
Janeiro - RJ, 22290-160
E-mail: mehecobrasil@
gmail.com

Indian Showroom

1460 M.I.E, Bahadurgarh,
Haryana-124507, INDIA
TEL: +91 96677 70747
E-mail: ind@saikangmedical.
com

China Headquarter

No.35 Lehong Road,
Zhangjiagang City,
Jiangsu Province, China.
TEL: +86 18021231901
E-mail: export@
saikangmedical.com

South East Asia Branch

**SAIKANG Medical
(Thailand) Co., Ltd.**
Ample Tower 909 Room 11/3
Debaratana Road,
Bangna Nuea, Bangkok
10260 Thailand
TEL: +66-936-569469

US Branch

**Healthward Technology,
INC**
17800 castleton st ste 180,
Los Angeles
TEL: +1 (380) 254 4913

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue